

Ibuprofen on LUX Amylose-1 5u in RP


Column: Lux® 5 µm Amylose-1, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4732-E0
Elution Type: Isocratic
Eluent A: 0.1% Acetic Acid in Water
Eluent B: 0.1% Acetic Acid in Acetonitrile
Gradient Profile:

Step No.	Time (min)	Pct A	Pct B
1	0	50	50

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)



Products used in this application:

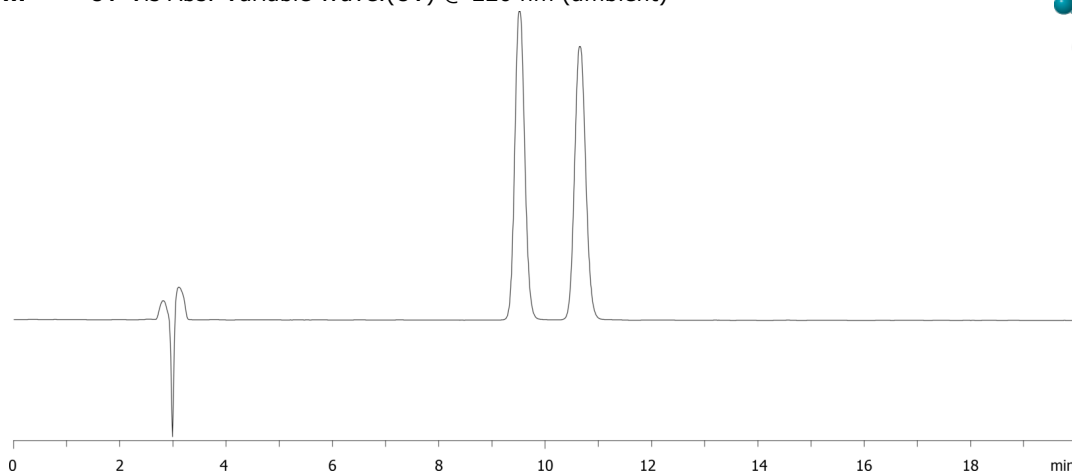

ANALYTES:

- 1 Ibuprofen
Retention Time: 9.495 min
- 2 Ibuprofen
Retention Time: 10.628 min

