


Ethanolamines using ZB-5MSplus 30 x 0.25 x 1.00**Column:** Zebron[™] ZB-5MSplus, GC Cap. Column 30 m x 0.25 mm x 1.00 μ m, Ea**Phase:****Dimensions:** 30 meters x 0.25 mm x 1 μ m**Order No:** 7HG-G030-22**Oven Profile:** 30 °C to 300 °C at 40 °C/min**Carrier Gas:** Constant Flow Helium, 1.4 mL/min**Injection:** Split 200:1 1 μ L @ 250°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Sim mode: Monoethanolamine (61,42), Diethanolamine (74, 56), Triethyleneglycol monomethylether (45,59,84),**Products used in this application:****ANALYTES:**

- 1 Monoethanolamine
- 2 Diethanolamine
- 3 Triethylene glycol monomethyl ether (I.S.)
- 4 Triethanolamine

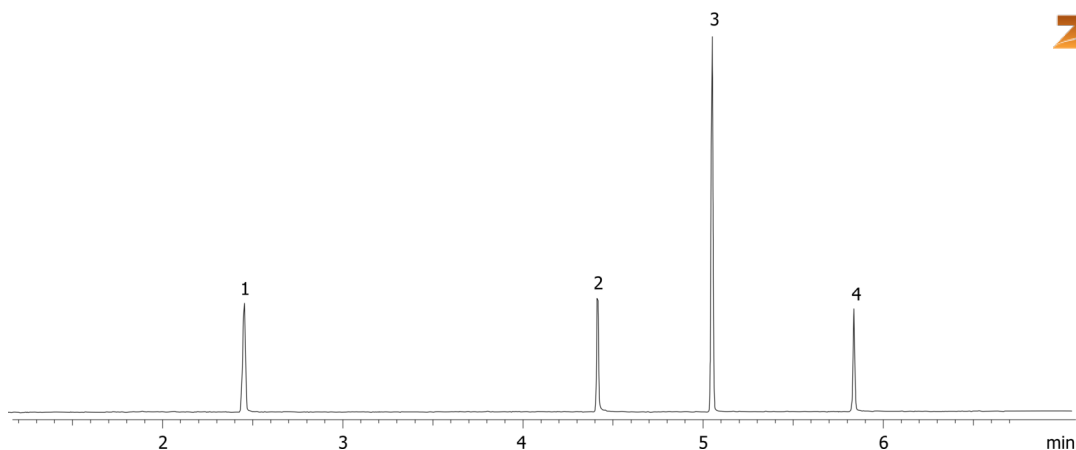