


Yarra 1.8 µm SEC-X150, Mouse Monoclonal Antibody

Column: N/A
Dimensions: N/A
Order No: N/A
Elution Type: N/A
Eluent A: 0.1 M Sodium Phosphate, pH 6.8
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 0.35 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

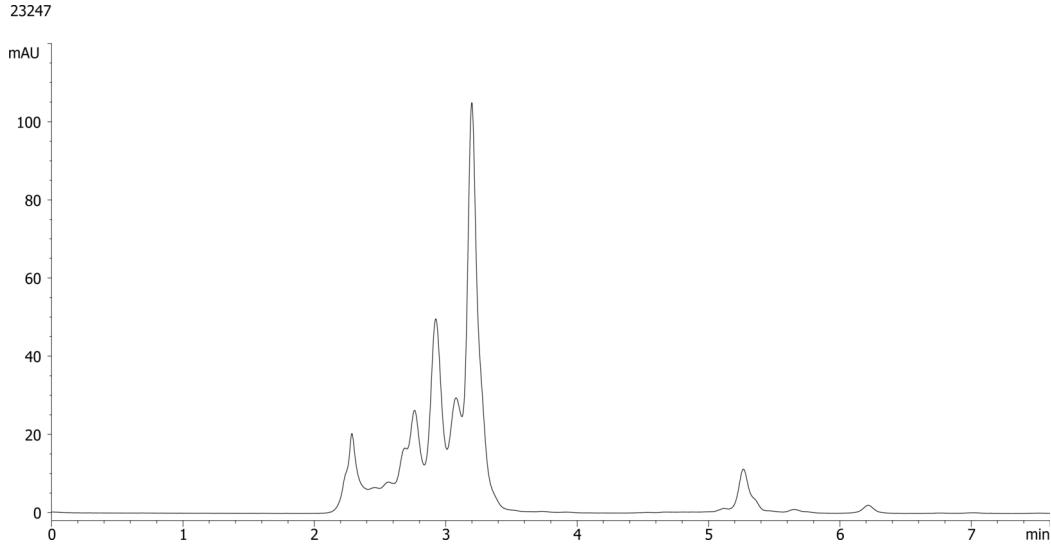

ANALYTES:

- 1 Mouse Monoclonal Antibody
Retention Time: 3.199 min

