

Fenoprofen on Lux 3u Amylose-1 250x4.6mm in SFC (MeOH and TFA)

Column: Lux® 5 µm Amylose-1, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4732-E0

Elution Type: Isocratic

Eluent A: CO₂


Eluent B: methanol with 0.5% trifluoroacetic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	2	80	20

Flow Rate: 3 mL/min

Col. Temp.: 30 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ (ambient)

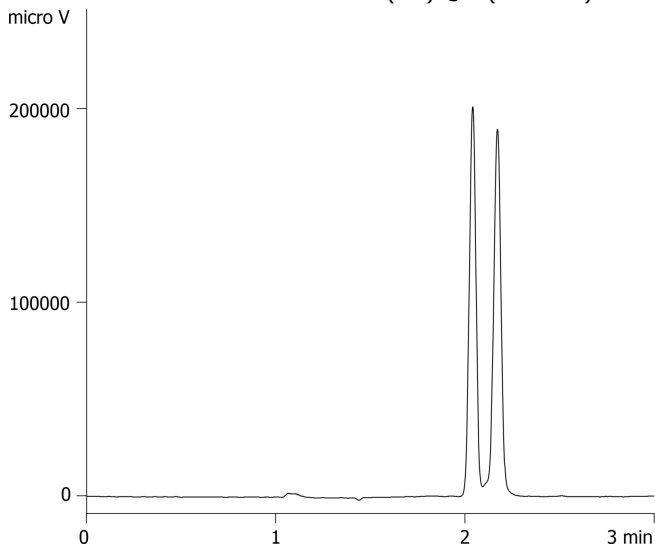

ANALYTES:



Products used in this application:

