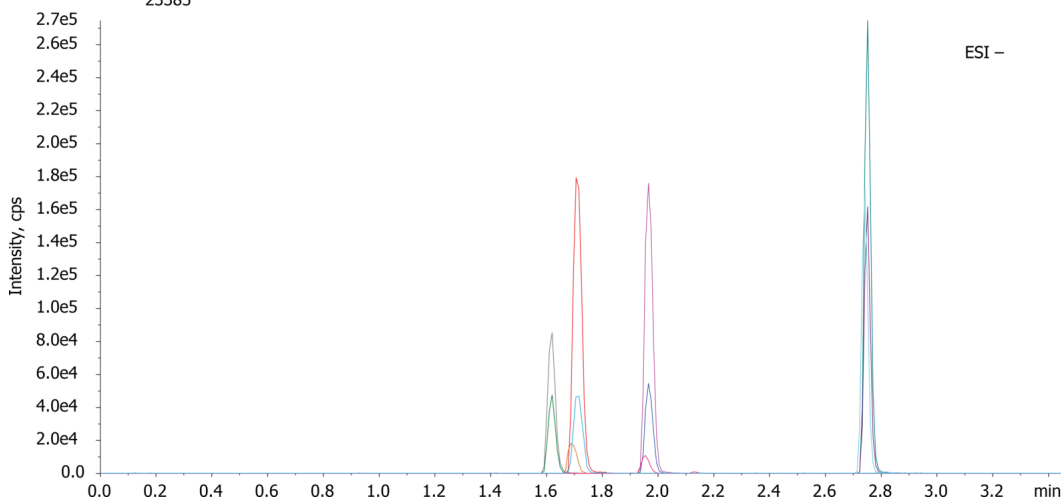


Fast Pain Panel ESI-**Column:** Kinetex® 2.6 µm Biphenyl 100 Å LC Column 30 x 2.1 mm, Ea**Dimensions:** 30 x 2.1 mm ID**Order No:** 00A-4622-AN**Elution Type:** Gradient**Eluent A:** 10mM ammonium formate in water**Eluent B:** 0.1% formic acid in methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	0.2	90	10
	3	2.5	10	90
	4	2.9	10	90
	5	3	90	10
	6	3.5	90	10

Flow Rate: 1 mL/min**Col. Temp.:** 35 °C**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)**Detector Info:** <a target="_blank"

href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

**ANALYTES:**

- 1 Butalbital
- 2 Phenobarbital
- 3 Secobarbital
- 4 THC-COOH
- 5 THC-COOH-D9
- 6 Secobarbital-D5
- 7 Butalbital-D5

