

Representative chromatogram for 5 g mussel extract

Column: Phenogel™ 5 µm 100A, LC Column 300 x 21.2 mm, Ea

Dimensions: 300 x 21.2 mm ID

Order No: 00H-0442-P0

Elution Type: Isocratic


Eluent A: Methylene Chloride

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

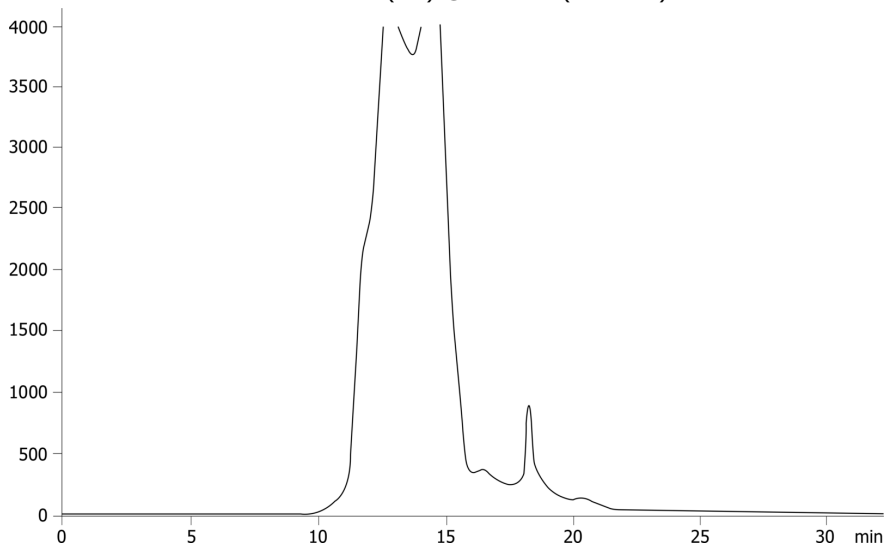

ANALYTES:

1 Mussel Extract



Products used in this application:

