

Nucleotides on Luna Omega Polar C18 50x2.1mm

Column: Luna Omega 1.6 um Polar C18 100 A LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4748-AN

Elution Type: Isocratic

Eluent A: 20mM Ammonium Formate pH 3.0

Gradient	Step No.	Time (min)	Pct A
Profile:	1	9	100


Flow Rate: 0.4 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 285 nm (ambient)



Products used in this application:

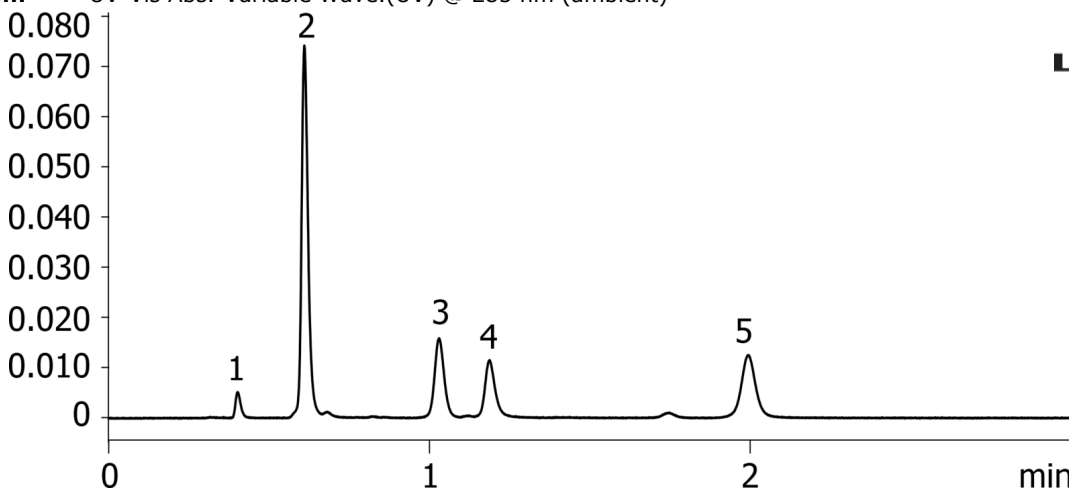

ANALYTES:

- 1 Thiourea
- 2 5-Fluorocytosine
- 3 Guanosine-5-Monophosphate
- 4 Adenine
- 5 Thymine

