

SAMHSA workplace drug testing panel using Luna Omega Polar C18 50x2.1mm

Column: Luna Omega 1.6 μ m Polar C18 100 A LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4748-AN

Elution Type: Gradient

Eluent A: Water with 0.1% formic acid

Eluent B: Acetonitrile with 0.1% formic acid


Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	85	15
	2	2	0	100
	3	3	0	100
	4	3.1	85	15

Flow Rate: 400 μ L/min

Col. Temp.: 25 $^{\circ}$ C

Detection: Tandem Mass Spec (MS-MS) @ (450 $^{\circ}$ C)

Detector Info: <a target="_blank"

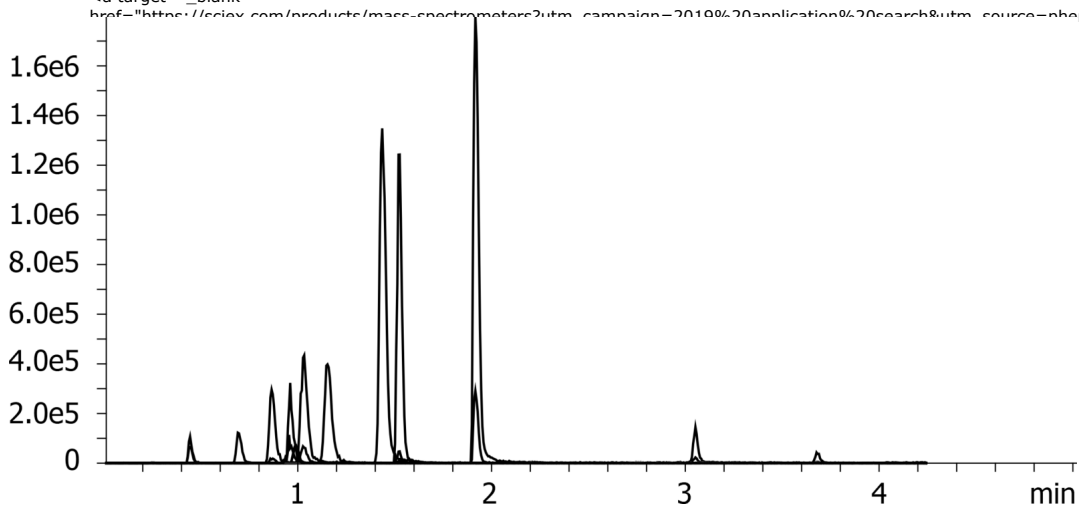

Products used in this application:



SAMHSA workplace drug testing panel using Luna Omega Polar C18 50x2.1mm

ANALYTES:

- 1** Morphine
- 2** Codeine
- 3** Amphetamine
- 4** MDA
- 5** Methamphetamine
- 6** 6-MAM
- 7** MDMA
- 8** MDEA
- 9** BZE
- 10** PCP
- 11** THC-COOH
- 12** THC

