


EPA 504 by GC/ECD on ZB-CLPesticides-2 (Hydrogen Gas)**Column:** Zebtron ZB-CLPesticides-2 Column 30m x 0.32mm x 0.25µm , Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G029-11**Oven Profile:** 30 °C for 2 min to 200 °C at 30 °C/min**Carrier Gas:** Constant Flow Hydrogen, 60 cm/sec**Injection:** Splitless :1 1 µL @ 250°C**Detection:** Electron Capture (ECD) (330°C)**Analyst Note:** Analytes are 1000 ng/mL in Hexane
Columns are split using a Y-splitter connected to a 5 m x 0.32 ID guard column (7AM-G000-00-GZ0)**Products used in this application:**

23635

**ANALYTES:**

- 1 Bromodichloromethane
- 2 Chlorodibromomethane
- 3 1,2-Dibromomethane (EDB)
- 4 1,1,1,2-Tetrachloroethane
- 5 Bromoform
- 6 1,2,3-Trichloropropane
- 7 1,2-Dibromo-3-chloropropane (DBCP)

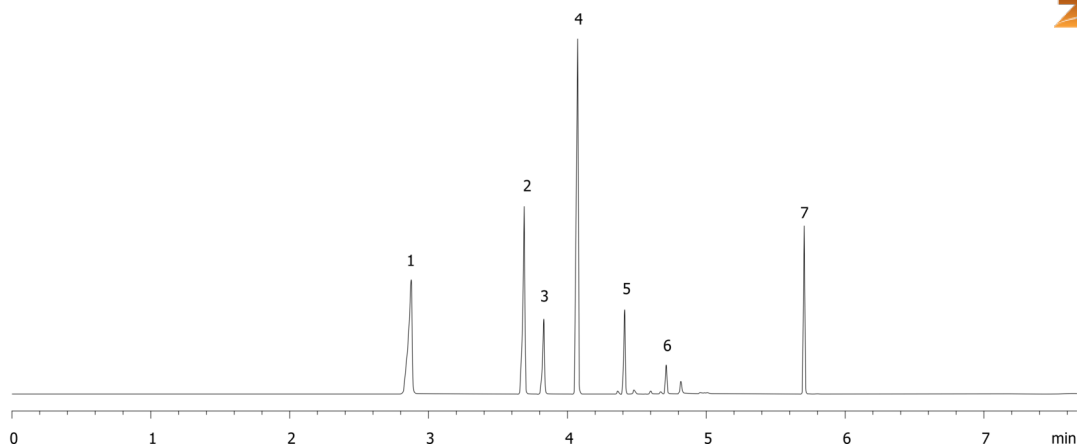