


OF tech note Fig 2, Negative ESI

Column: Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 3.0 mm, Ea**Dimensions:** 50 x 3 mm ID**Order No:** 00B-4622-Y0**Elution Type:** Gradient**Eluent A:** 10 mM ammonium formate in water**Eluent B:** Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	4	5	95
	3	5	5	95
	4	5.01	90	10
	5	7	90	10

Flow Rate: 500 µL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)**Detector Info:** <a target="_blank"

App ID 23720

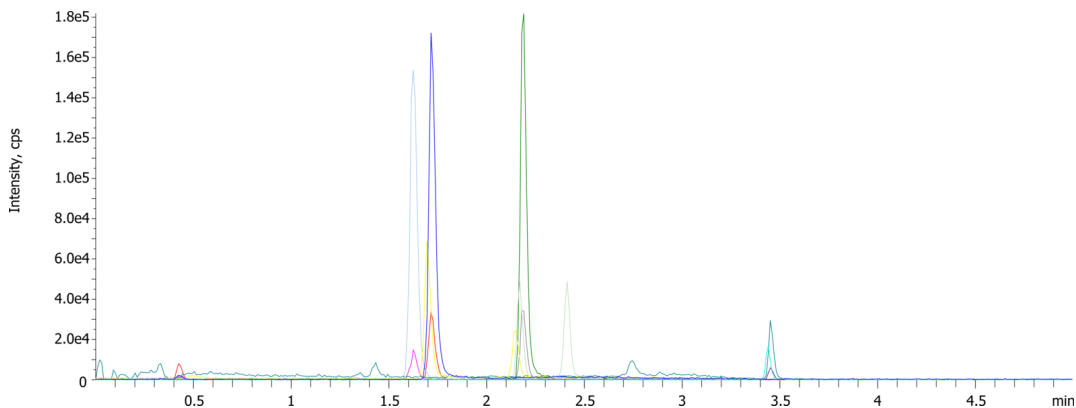

Products used in this application:



ANALYTES:

1 Pain Panel

