


Follitropin System Suitability on 5 µm Hydrophilic Silica

Column: BioSep™ 5 µm SEC-s2000 145 Å, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-2145-K0
Elution Type: Isocratic
Eluent A: 100 mM Sodium Phosphate, 200 mM Sodium Sulfate, pH 6.7
Gradient


Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 0.5 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 215 nm (ambient)

23789

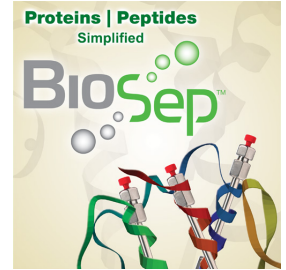

Products used in this application:

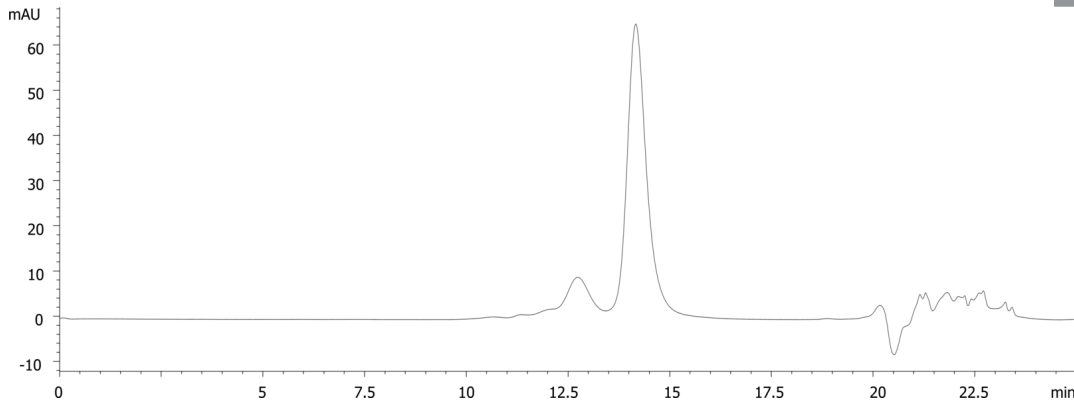

ANALYTES:

- 1 BSA
Retention Time: 12.736 min
- 2 Follitropin
Retention Time: 14.161 min

