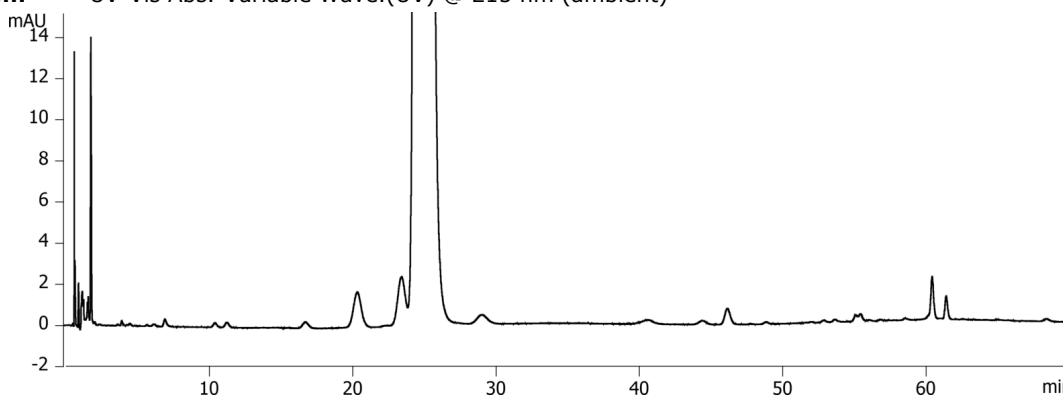


Ezetimibe USP on Luna Omega 5u Polar C18 150x4.6mm**Column:** Luna Omega 5u Polar C18 100A, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4754-E0**Elution Type:** Gradient**Eluent A:** Water:Methanol (90:10)**Eluent B:** Acetonitrile:Metanol (90:10)

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	68	32
	2	37	68	32
	3	60	45	55
	4	70	45	55
	5	80	15	85
	6	90	15	85
	7	90.1	68	32
	8	100	68	32

Flow Rate: 2 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 215 nm (ambient)**ANALYTES:****1** Ezetembe**Products used in this application:**