

FAMES In Regular Olive Oil by GC/FID on Zebron ZB-FAME (30 m)

Column: Zebron ZB-FAME, GC Cap. Column 30m x 0.25mm x 0.2um , Ea
Phase: Proprietary Pesticides Phase
Dimensions: 30 meters x 0.25 mm x 0.2 µm
Order No: 7HG-G033-10
Oven Profile: 100 °C for 2 min to 140 °C @ 10 °C/min to 190 °C @ 3 °C/min to 260 °C @ 30 °C/min for 2 min
Carrier Gas: Constant Flow Helium, 1.2 mL/min
Injection: Split 50:1 1 µL @ 240°C
Detection: Flame Ionization (FID) (260°C)
Analyst Note: Recommended Liner: Zebron PLUS Single Taper with Wool, 4 mm ID
 Liner Part No.: AG2-0A11-05 (for Agilent systems)
 Inlet Seal: AG0-8620 (Gold Plated Easy Seal)
 Septum: AG0-4696 (PhenoRed-400)




Products used in this application:



Sample Preparation:

1. Strata® Si-1 Tube, 1 g/6 mL (Part No.: 8B-S012-JCH) on a vacuum or positive pressure manifold
2. Wash cartridge with 6 mL of hexane
3. Load oil solution (0.12 g of oil in 0.5 mL of hexane)
4. Elute with 10 mL of hexane/diethyl ether (87:13)
5. Evaporate eluate under a steady stream of nitrogen
6. Dissolve purified oil residue in 1mL of heptane
7. Add 0.1 mL of 2N potassium hydroxide in methanol to purified oil solution
8. Cap tube and shake vigorously for 15 seconds
9. Leave to separate until upper layer becomes clear
10. Extract upper layer for GC analysis

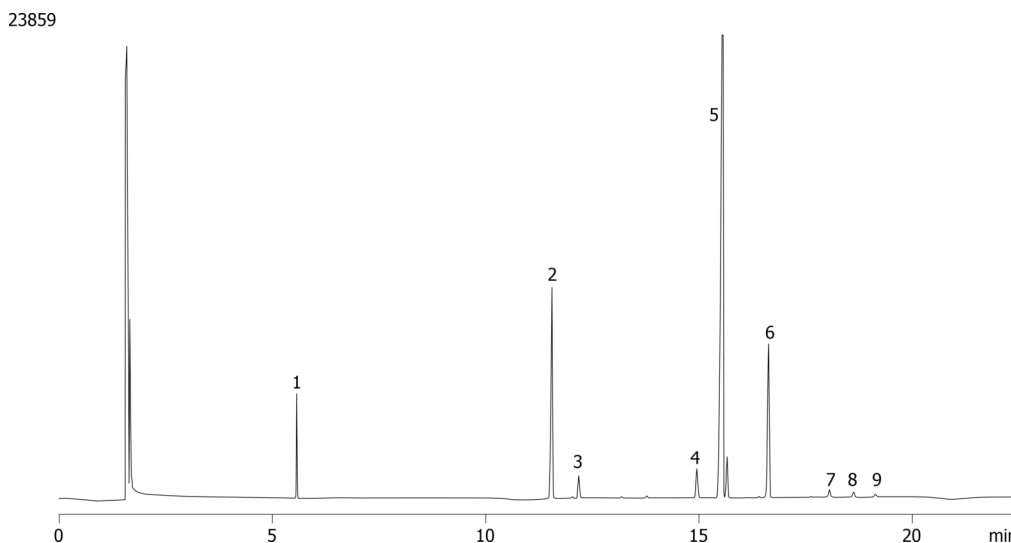

ANALYTES:

- 1 C11:0 (I.S.)
- 2 C16:0
- 3 C16:1 cis 9
- 4 C18:0
- 5 C18:1 cis 9
- 6 C18:2 cis 9,12
- 7 C18:3 cis 9,12,15
- 8 C20:0
- 9 C20:1 cis 11

©2005 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com



Phenomenex products are available worldwide.

www.phenomenex.com.cn

info@phenomenex.com

Sample Preparation Details

for GC

Application ID No.: 23859

FAMES In Regular Olive Oil by GC/FID on Zebron ZB-FAME (30 m)

PRODUCT DESCRIPTION:

Strata® SI-1 Silica (55 µm, 70 Å), 1 g / 6 mL, Tubes , 30/Pk

Order No.: 8B-S012-JCH

SOLID PHASE EXTRACTION (SPE) PROCEDURE:

Note: The solvent volumes shown below are for a 1 g bed mass.

The solvent volumes will need to be adjusted for a smaller or larger bed mass.

Condition:

Load:

Add 0.12 g of oil to 0.5 mL of Hexane, load onto cartridge

Wash:

Dry:

Elute:

Final Prep and Analysis:

To reconstituted sample, add 0.1 mL of 2 N Potassium hydroxide in methanol. Cap tube and shake vigorously for 15 seconds. Leave to separate until upper layer becomes clear. Extract

Inject: 1 µL on HPLC Flame Ionization (FID) @ (260°C)

ANALYTES:	Spiked Conc. (ng/mL)	Log P	pKa	% Rec	%RSC (n=0)
1 C11:0 (I.S.)	0				
2 C16:0	0				
3 C16:1 cis 9	0				
4 C18:0	0				
5 C18:1 cis 9	0				
6 C18:2 cis 9,12	0				
7 C18:3 cis 9,12,15	0				
8 C20:0	0				
9 C20:1 cis 11	0				

Note: This method is designed as a convenient starting point for further investigation and can be tailored to meet your extraction goals. Call your local Phenomenex Representative for assistance in method development and optimization techniques.

©2025 Phenomenex Inc. All rights reserved.

For more information contact your Phenomenex Representative at info@phenomenex.com



Phenomenex products are available worldwide.

www.phenomenex.com.cn

info@phenomenex.com