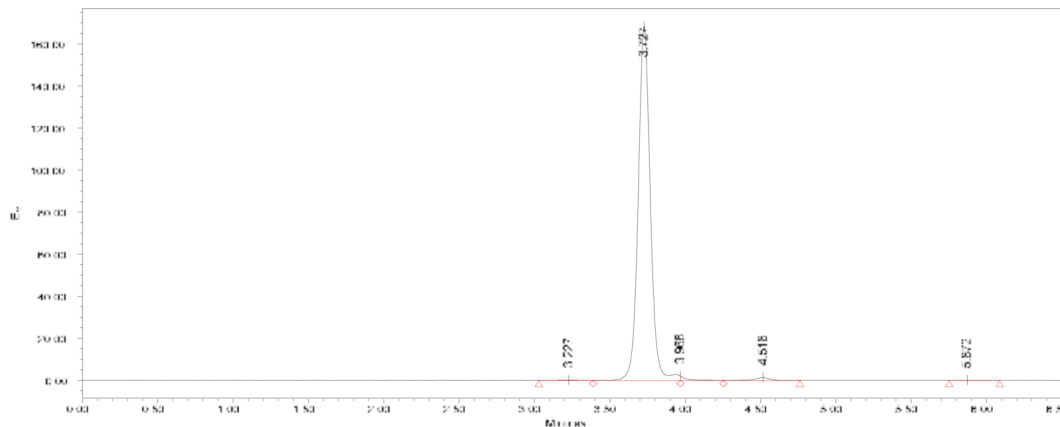


Adalimumab on Yarra SEC-X300**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 50 mM Potassium Phosphate, 250 mM KCl, pH 6.8

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.35 mL/min**Col. Temp.:** 25 °C**Detection:** Fluorescence (FLUOR) @ 280 nm (25 °C)

23903

**ANALYTES:****1** Adalimumab