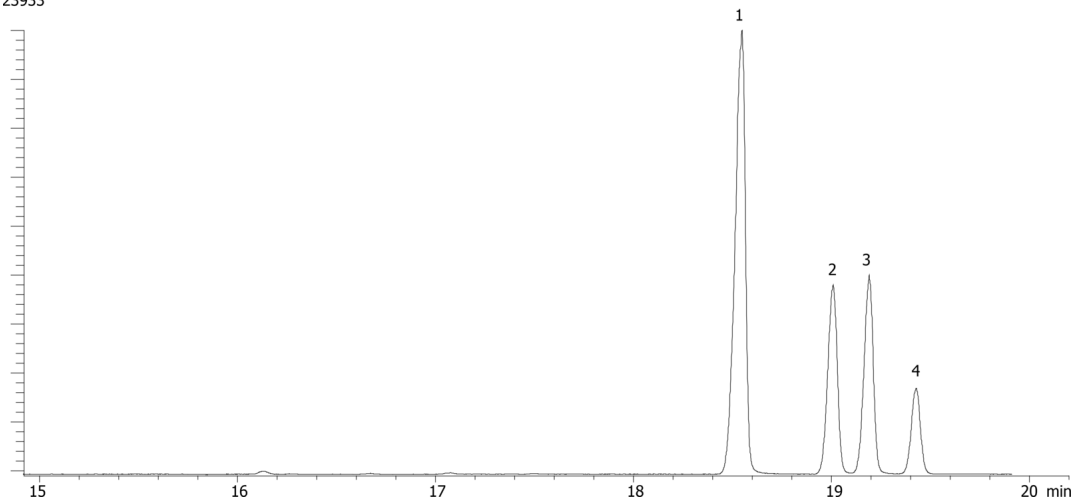


Vit B12 - Adenosylcobalamin Product Ion Tuning File**Column:** Synergi™ 2.5 µm Fusion-RP 100 Å, LC Column 50 x 2 mm, Ea**Dimensions:** 50 x 2 mm ID**Order No:** 00B-4423-B0**Elution Type:** Gradient**Eluent A:** 10 mM Ammonium Formate pH 3.5**Eluent B:** MeOH

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	98	2
	2	0.2	98	2
	3	1.8	2	98
	4	4	2	98
	5	4.01	98	2
	6	5.5	98	2

Flow Rate: 350 µL/min**Col. Temp.:** 21 °C**Detection:** Tandem Mass Spec (MS-MS) @ (600 °C)**Detector Info:** SCIEX<

23933

**ANALYTES:****1** Adenosylcobalamin

Retention Time: 2.5 min

