

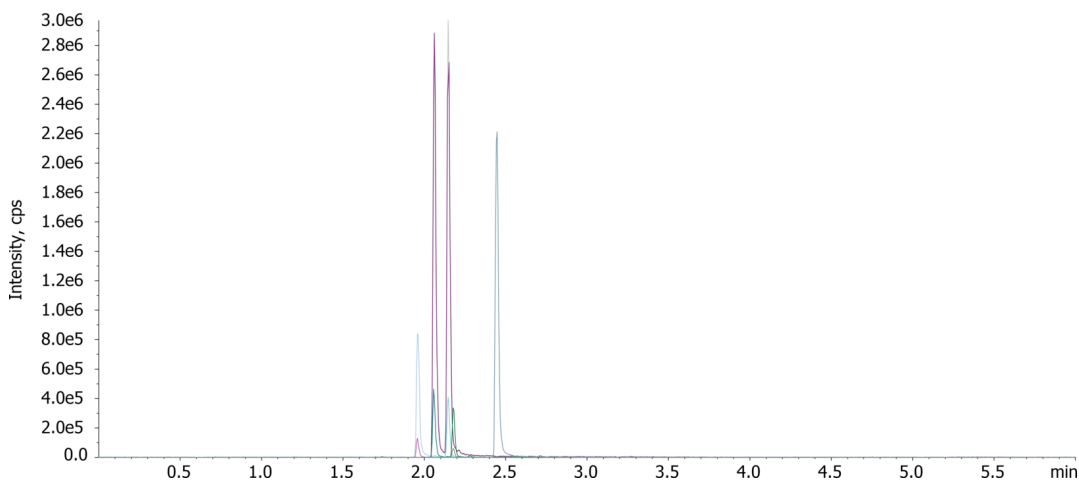
Bath salts (250ng/ml) from urine using Kinetex 2.6u C18 50x4.6mm - Extracted sample**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 4.6 mm, Ea**Dimensions:** 50 x 4.6 mm ID**Order No:** 00B-4462-E0**Elution Type:** Gradient**Eluent A:** Water with 0.1% formic acid**Eluent B:** AcCN with 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	4	5	95
	3	4.1	95	5
	4	6	95	5

Flow Rate: 700 mL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)**Detector Info:** <a target="_blank"

href="https://www.phenomenex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

23941

**ANALYTES:**

- 1 Methylone
Retention Time: 1.96 min
- 2 Ethylone
Retention Time: 2.06 min
- 3 Butylone
Retention Time: 2.15 min
- 4 Mephedrone
Retention Time: 2.18 min
- 5 Alpha-PVP
Retention Time: 2.44 min

