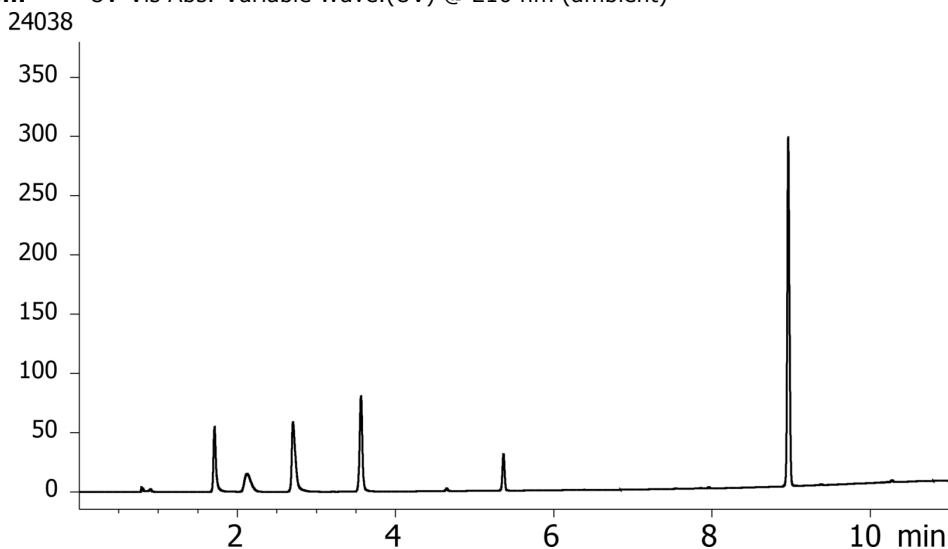


Water Soluble Vitamins on Kinetex 2.6u Polar C18 100x4.6mm**Column:** Kinetex® 2.6 µm Polar C18 100 Å, LC Column 100 x 4.6 mm, Ea**Dimensions:** 100 x 4.6 mm ID**Order No:** 00D-4759-E0**Elution Type:** Gradient**Eluent A:** 20mM Potassium Phosphate pH 3.5**Eluent B:** Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	1	100	0
	3	10	40	60
	4	11	40	60
	5	11.1	100	0
	6	15	100	0

Flow Rate: 1.2 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)**ANALYTES:**

- 1 Thiamine
- 2 Nicotinamide
- 3 Pyridoxal
- 4 Pyridoxine
- 5 Panthothenic Acid
- 6 Riboflavin

**Products used in this application:**