

EPA 8330A on Kinetex 2.6u Polar C18 100x4.6mm

Column: Kinetex® 2.6 µm Polar C18 100 Å, LC Column 100 x 4.6 mm, Ea

Dimensions: 100 x 4.6 mm ID

Order No: 00D-4759-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Methanol

Gradient


Step No.	Time (min)	Pct A	Pct B
1	12	72	28

Profile:

Flow Rate: 1.35 mL/min

Col. Temp.: 40 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

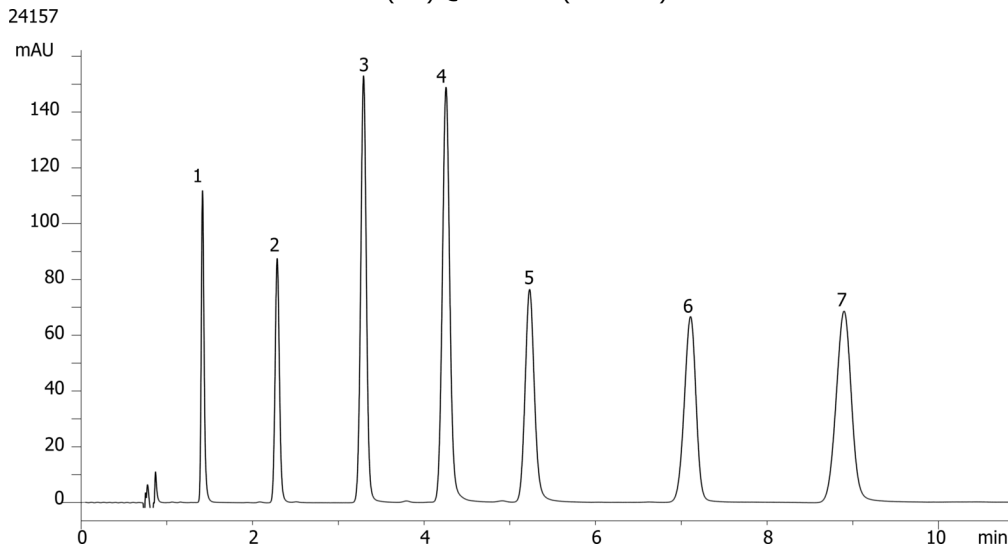

ANALYTES:

- 1 HMX
- 2 RDX
- 3 1,3,5-Trinitrobenzene
- 4 1,3-Dinitrobenzene
- 5 Nitrobenzene
- 6 2,4,6-Trinitrotoluene
- 7 2,4-Dinitrotoluene



Products used in this application:

