

12 Cannabinoids in AcCN using Kinetex 2.6µm Polar C18 150x4.6mm

Column: Kinetex® 2.6 µm Polar C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4759-E0

Elution Type: Gradient

Eluent A: Water with 0.1% phosphoric acid

Eluent B: Acetonitrile with 0.1% phosphoric acid

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	25	75
	2	6	0	100


Flow Rate: 1.8 mL/min

Col. Temp.: 50 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)



Products used in this application:

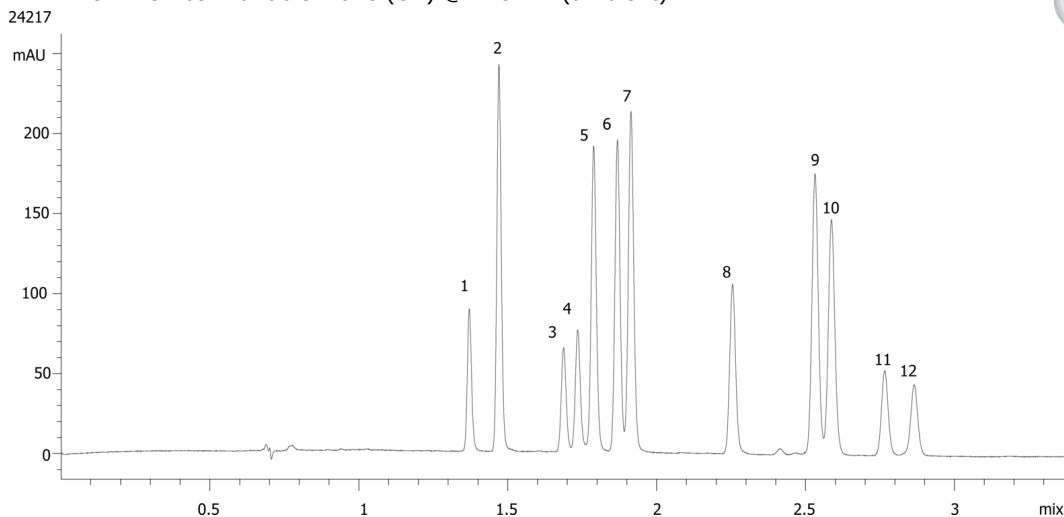

12 Cannabinoids in AcCN using Kinetex 2.6um Polar C18 150x4.6mm

ANALYTES:

- 1** CBDV-A
- 2** CBDV
- 3** CBD-A
- 4** CBG-A
- 5** CBG
- 6** CBD
- 7** THC-V
- 8** CBN
- 9** Delta-9-THC
- 10** Delta-8-THC
- 11** CBC
- 12** THC-A

