

Organic Acids on Kinetex 2.6u Polar C18 150x4.6mm

Column: Kinetex® 2.6 µm Polar C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4759-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium Phosphate pH 1.5

Gradient

Step No.	Time (min)	Pct A
1	0	100


Profile:

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

App Id 24268

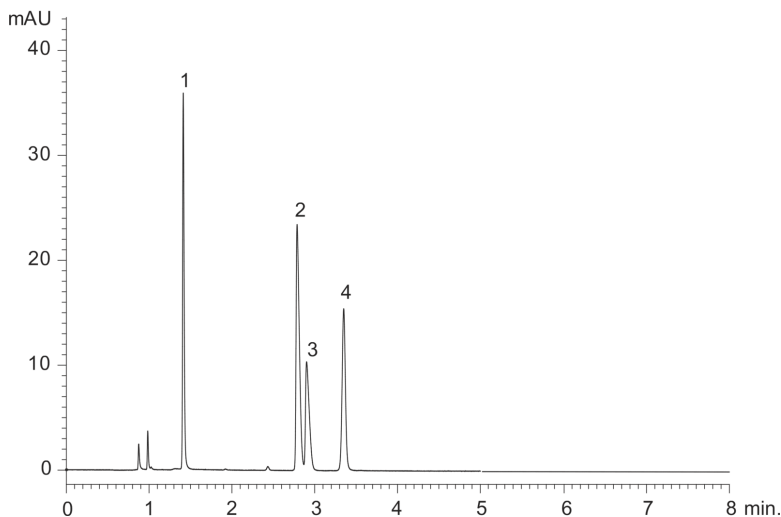

ANALYTES:

- 1 Malic Acid
- 2 Citric Acid
- 3 Fumaric Acid
- 4 Succinic Acid



Products used in this application:

