

## Epinephrine on Lux 5 $\mu$ Amylose-1 250x4.6mm

**Column:** Lux® 5  $\mu$ m Amylose-1, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4732-E0

**Elution Type:** Isocratic

**Eluent A:** n-Hexane + 0.2% TFA + 0.2% DEA

**Eluent B:** 2-Propanol + 0.2% TFA + 0.2% DEA

**Eluent C:** Methanol + 0.2% TFA + 0.2% DEA

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B | Pct C |
|-------------------|----------|------------|-------|-------|-------|
| 1                 | 0        | 90         | 5     | 5     |       |


**Flow Rate:** 1 mL/min

**Col. Temp.:** 30 °C

**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 280 nm (ambient)



Products used in this application:

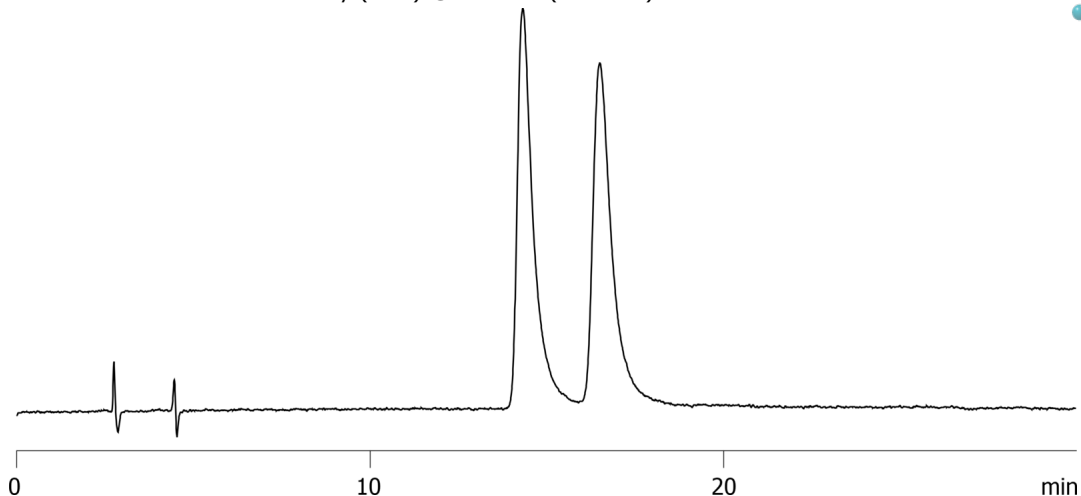

### ANALYTES:

1 ( $\Delta\pm$ )-Epinephrine Hydrochloride

