


Chiral Separation of Nadolol on Lux 5 μ m i-Amylose-1 in NP**Column:** Lux[®] 5 μ m i-Amylose-1, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4762-E0**Elution Type:** Isocratic**Eluent A:** 0.1% DEA in Hexane**Eluent B:** 0.1% DEA in IPA

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	85	15

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)**Products used in this application:****ANALYTES:**

- 1 Nadolol
Retention Time: 12.355 min
- 2 Nadolol
Retention Time: 12.915 min
- 3 Nadolol
Retention Time: 22.312 min
- 4 Nadolol
Retention Time: 26.359 min

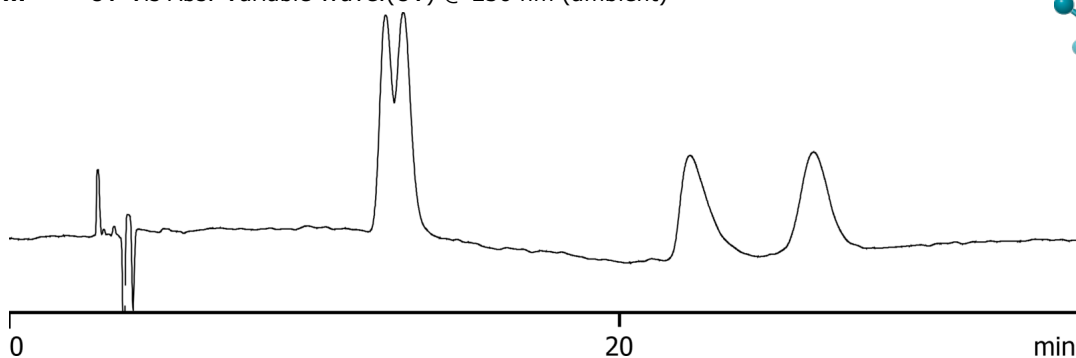