

Chiral Separation of 2-Bromo-1-Indanone on Lux 5µm i-Amylose-1 in Dichloromethane


Column: Lux® 5 µm i-Amylose-1, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4762-E0
Elution Type: Isocratic
Eluent A: 0.1% TFA in Hexane
Eluent B: 0.1% TFA in Dichloromethane
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	65	35

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)



Products used in this application:

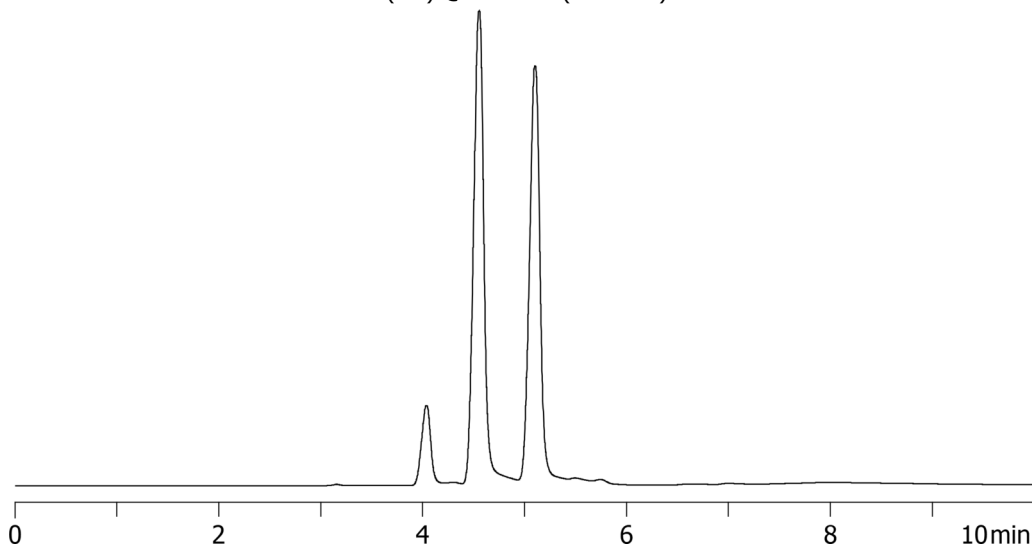

ANALYTES:

- 1 2-Bromo-1-Indanone
Retention Time: 4.552 min
- 2 2-Bromo-1-Indanone
Retention Time: 5.102 min

