

Cannabinoids (12) using Kinetex 2.6µm Polar C18 150x4.6mm

Column: Kinetex® 2.6 µm Polar C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4759-E0

Elution Type: Gradient

Eluent A: Water with 0.1% TFA

Eluent B: Acetonitrile with 0.1% TFA


Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	25	75
	2	7	0	100

Flow Rate: 1 mL/min

Col. Temp.: 55 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)

App ID 24347

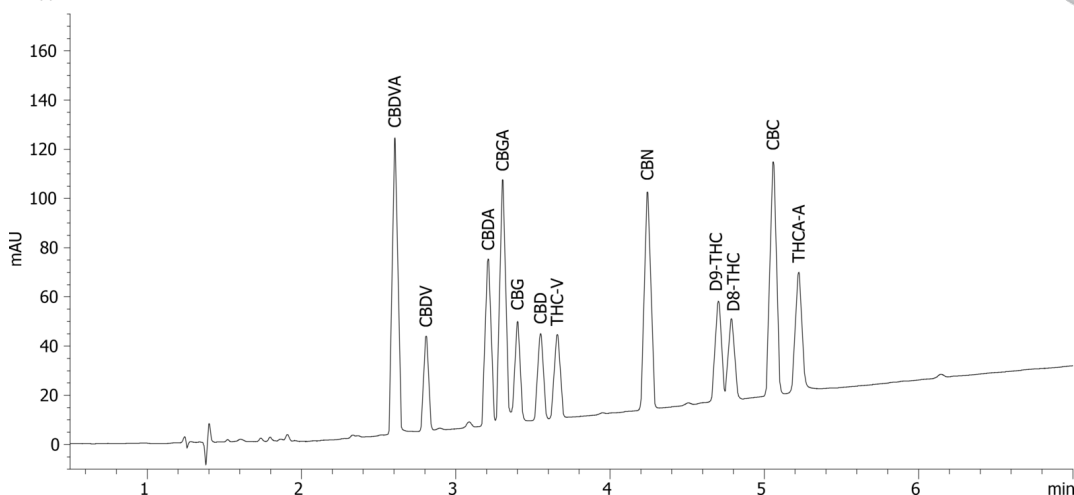

Products used in this application:



Cannabinoids (12) using Kinetex 2.6um Polar C18 150x4.6mm

ANALYTES:

- 1** CBDVA
- 2** CBDV
- 3** CBDA
- 4** CBGA
- 5** CBG
- 6** CBD
- 7** THC-V
- 8** CBN
- 9** D9-THC
- 10** D8-THC
- 11** CBC
- 12** THCA-A

