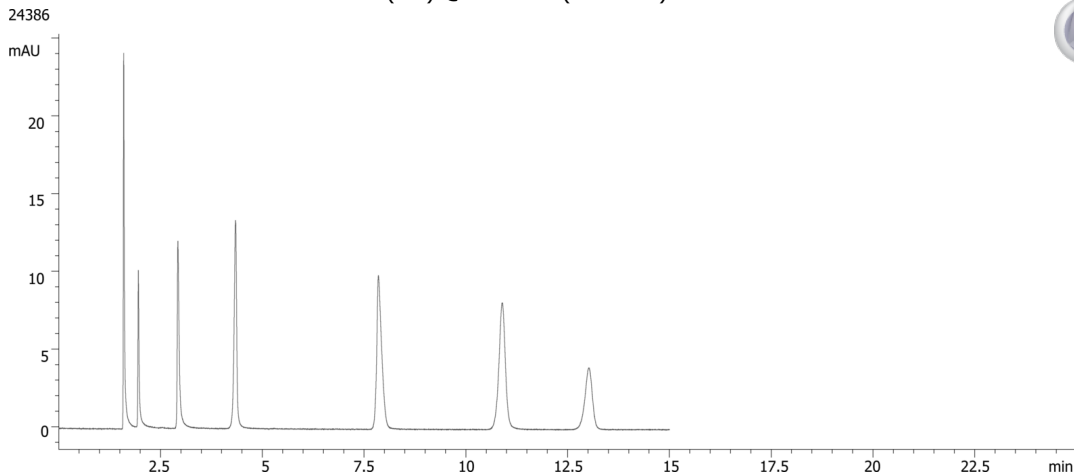


Kinetex C18 Selectivity Tech Note - Fig 3 Kinetex XB-C18**Column:** Kinetex® 2.6 µm XB-C18 100 Å, LC Column 50 x 4.6 mm, Ea**Dimensions:** 50 x 4.6 mm ID**Order No:** 00B-4496-E0**Elution Type:** Isocratic**Eluent A:** 20mM Ammonium Acetate pH 6.9**Eluent B:** Methanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	98	2

Flow Rate: 1 mL/min**Col. Temp.:** 25 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)**Products used in this application:****ANALYTES:****1** Same as 24384