

Butter Matrix - Phenolic Antioxidants Using Kinetex 2.6u Polar C18

Column: Kinetex® 2.6 µm Polar C18 100 Å, LC Column 100 x 3.0 mm, Ea

Dimensions: 100 x 3 mm ID

Order No: 00D-4759-Y0

Elution Type: Gradient

Eluent A: Water + 0.5% Acetic Acid

Eluent B: Acetonitrile + 0.5% Acetic Acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25
	2	5	0	100
	3	15	0	100

Flow Rate: 0.4 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

24458



Products used in this application:

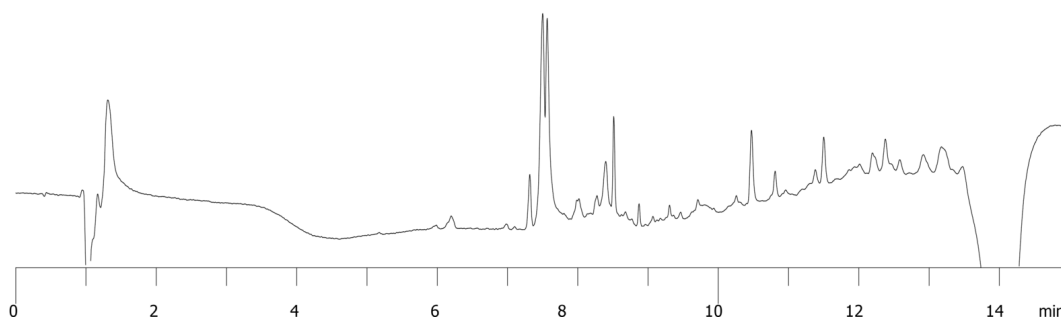

Butter Matrix - Phenolic Antioxidants Using Kinetex 2.6u Polar C18

ANALYTES:

- 1** Propyl Gallate
- 2** Ethoxyquin
- 3** TBHQ
- 4** NDGA
- 5** BHA
- 6** Ionox â€" 100
- 7** Vitamin E
- 8** Octyl Gallate
- 9** BHT
- 10** Lauryl Gallate
- 11** Vitamine E Acetate

