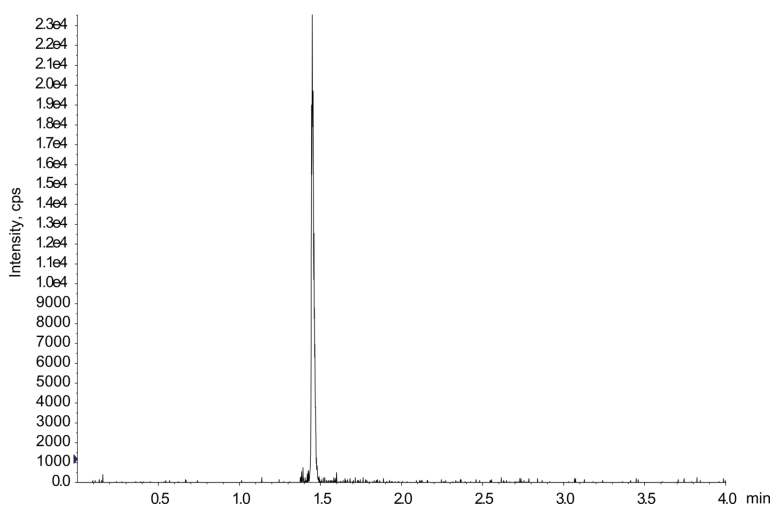


Naltrexone on Beta Gane with Kura Enzyme**Column:** Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4622-AN**Elution Type:** Gradient**Eluent A:** 0.1% Formic Acid in Water**Eluent B:** 01.% Formic acid in methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	2	0	100
	3	2.1	95	5
	4	4	95	5

Flow Rate: 0.5**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)

App ID 24462

**Products used in this application:****ANALYTES:****1** naltrexone

Retention Time: 1.45 min

