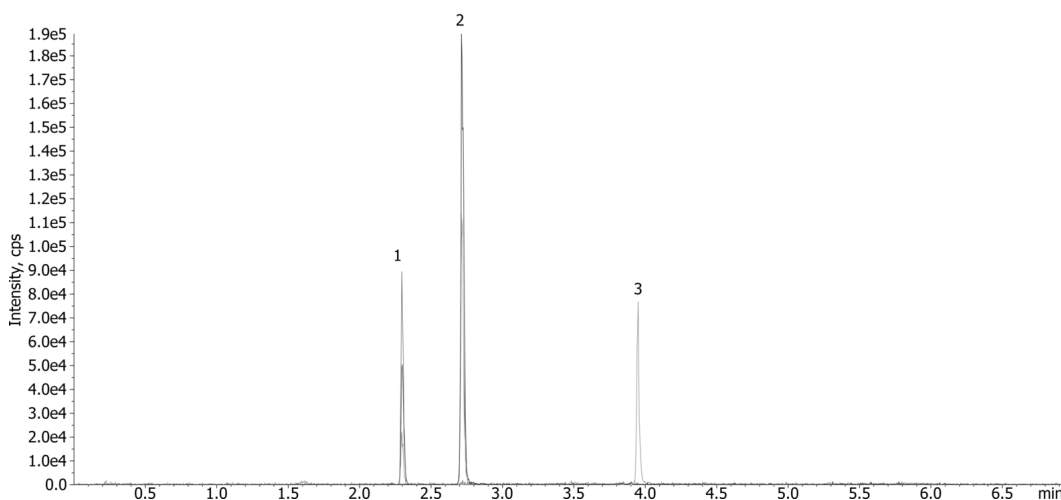


Nor/Buep + THC-COOH Novum ESI+**Column:** Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4622-AN**Elution Type:** Gradient**Eluent A:** 0.1% formic acid in water**Eluent B:** 0.1% formic acid in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	0	100
	3	5.1	95	5
	4	7	95	5

Flow Rate: 0.5 mL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)**Detector Info:** SCIEX<

24485

**ANALYTES:****1** Norbuprenorphine

Retention Time: 2.31 min

2 Buprenorphine

Retention Time: 2.71 min

3 THC-COOH

Retention Time: 3.95 min

**Products used in this application:**