

Chiral Separation of Suprofen on Lux 5µm i-Amylose-1 in Ethyl Acetate

Column: Lux® 5 µm i-Amylose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4762-E0

Elution Type: Isocratic

Eluent A: 0.1% DEA in Hexane

Eluent B: 0.1% TFA in Ethyl acetate

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	70	30	


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)
24488



Products used in this application:

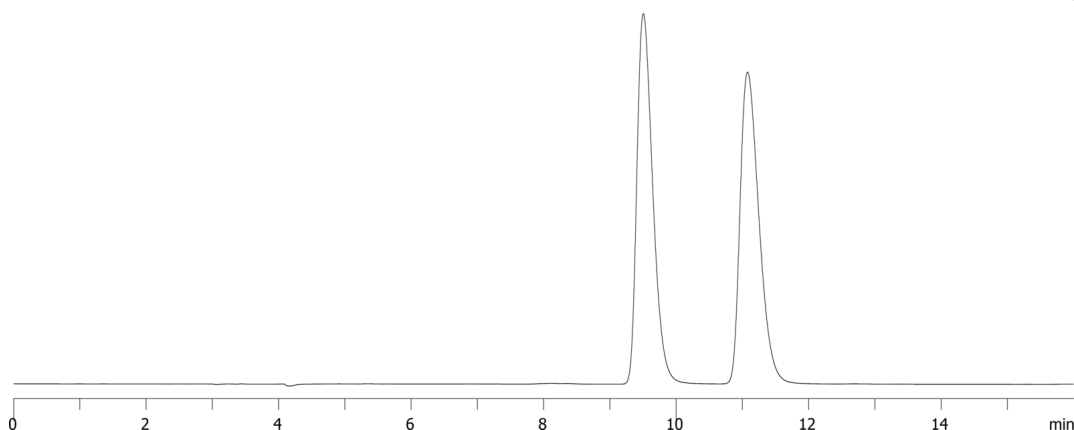

ANALYTES:

- 1 Suprofen
Retention Time: 9.503 min
- 2 Suprofen
Retention Time: 11.076 min

