

Chiral Separation of Chlorthalidone on Lux 5µm i-Amylose-1 in NP

Column: Lux® 5 µm i-Amylose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4762-E0

Elution Type: Isocratic

Eluent A: MTBE

Eluent B: Ethanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	1	0	90	10


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)
24499



Products used in this application:

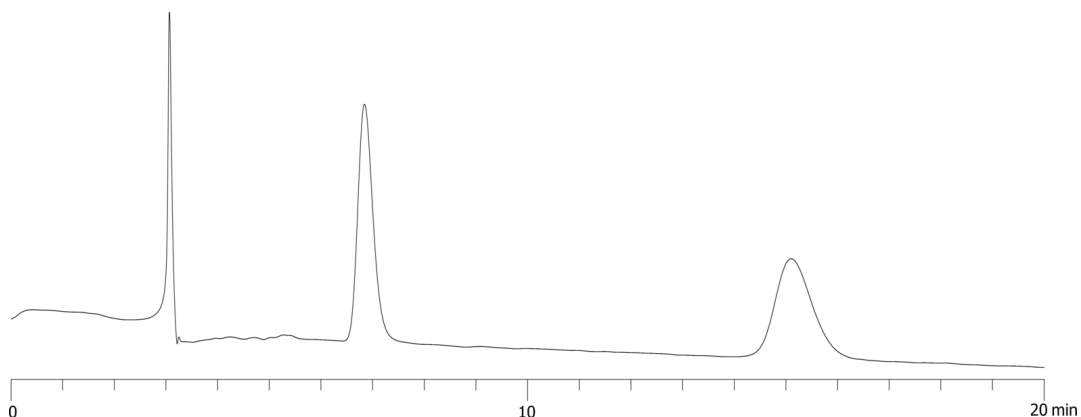

ANALYTES:

- 1 Chlorthalidone
Retention Time: 6.841 min
- 2 Chlorthalidone
Retention Time: 15.1 min

