


**Cold Pressed Orange Oil on ZB-1plus****Column:** Zebron ZB-1plus Capillary GC, Column 10m x 0.10mm IDx 0.1um, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 10 meters x 0.1 mm x 0 µm**Order No:** 7CB-G031-02**Oven Profile:** 60°C to 130°C at 10°C/min to 280°C at 30°C for 3min**Carrier Gas:** Constant Flow Helium, 0.3 mL/min**Injection:** Split 120:1 0.2 µL @ 160°C**Detection:** Mass Selective (MSD) (0°C)**Analyst Note:** Sample was 10% in dichloromethane**Products used in this application:**

## Cold Pressed Orange Oil on ZB-1plus

### ANALYTES:

- 1  $\alpha$ -Pinene
- 2 Sabinene
- 3 beta-Myrcene
- 4 Octanal
- 5 3-Carene
- 6 Limonene
- 7 Nonanal
- 8 Linalool
- 9 cis-Limonene oxide
- 10 trans-Limonene oxide
- 11 Citronellal
- 12  $\alpha$ -Terpineol
- 13 Decanal
- 14 Carvone
- 15 Neral
- 16 Geranial
- 17  $\alpha$ -Cubebene
- 18  $\beta$ -Cubenene
- 19 Dodecanal
- 20 Valencene
- 21 Cadinene
- 22 Nootkatone

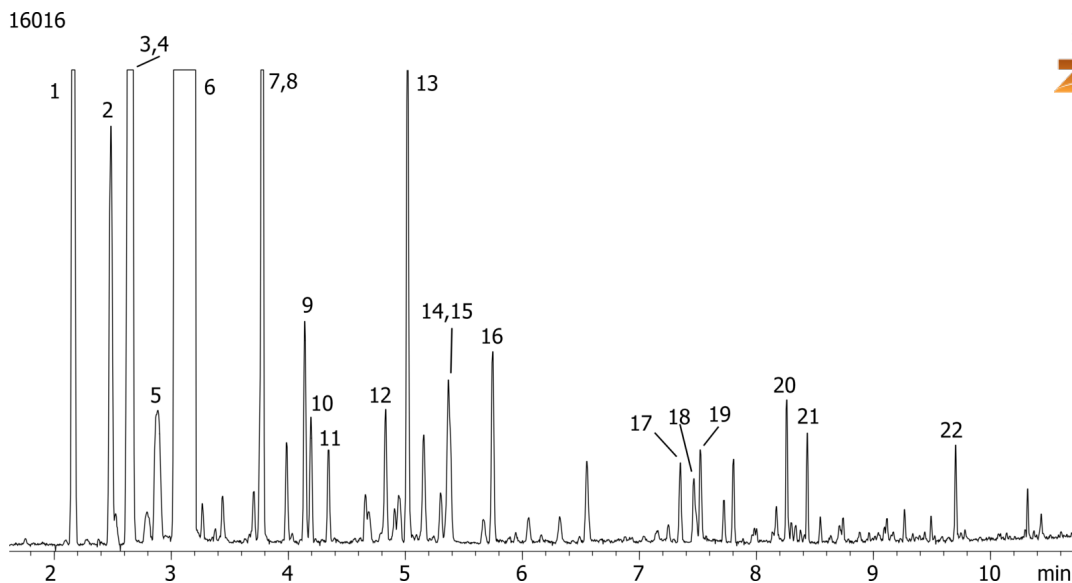