


Phenols ZB-5plus**Column:** Zebtron ZB-5plus Capillary GC, Column 30m x 0.25mm x 0.25um, Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G032-11**Oven Profile:** 60 °C to 140 °C @ 5 °C/min to 280 °C @ 10 °C/min.**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 5:1 1 µL @ 240°C**Detection:** Mass Selective (MSD) (230°C)**Analyst Note:** 50ppm in dichloromethane**Products used in this application:**

Phenols ZB-5plus

ANALYTES:

- 1 Phenol
- 2 2-Chlorophenol
- 3 2-Methylphenol
- 4 4-Methylphenol
- 5 3-Methylphenol
- 6 2,6-Dimethylphenol
- 7 2-Nitrophenol
- 8 2-Ethylphenol
- 9 2,4-Dimethylphenol
- 10 3,5-Dimethylphenol
- 11 2,5-Dimethylphenol
- 12 4-Ethylphenol
- 13 3-Ethylphenol
- 14 2,4-Dichlorophenol
- 15 Benzoic acid
- 16 2,3-Dimethylphenol
- 17 3,4-Dimethylphenol
- 18 2,6-Dichlorophenol
- 19 4-Chloro-3-methylphenol
- 20 2,4,6-Trichlorophenol
- 21 2,4,5-Trichlorophenol
- 22 2,4-Dinitrophenol
- 23 4-Nitrophenol
- 24 2,3,4,6-Tetrachlorophenol
- 25 4,6-Dinitro-2-methylphenol
- 26 Pentachlorophenol
- 27 Dinoseb

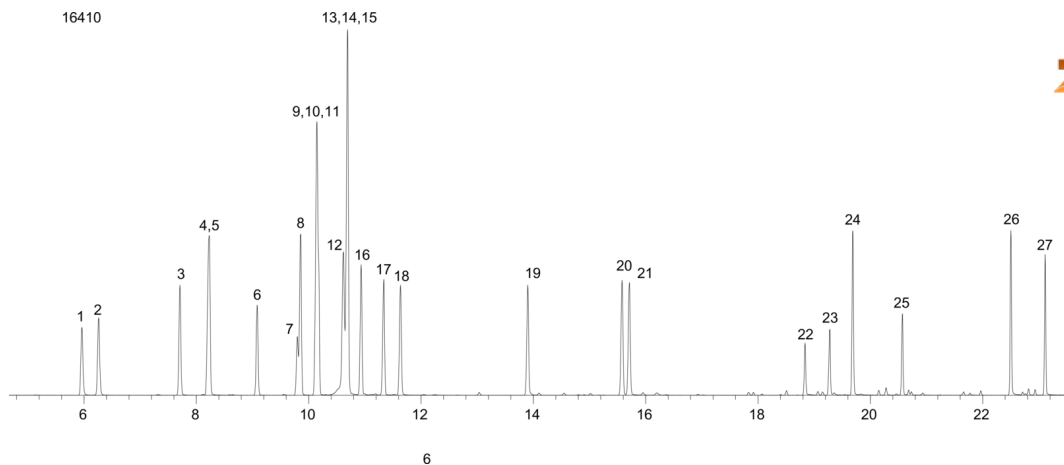