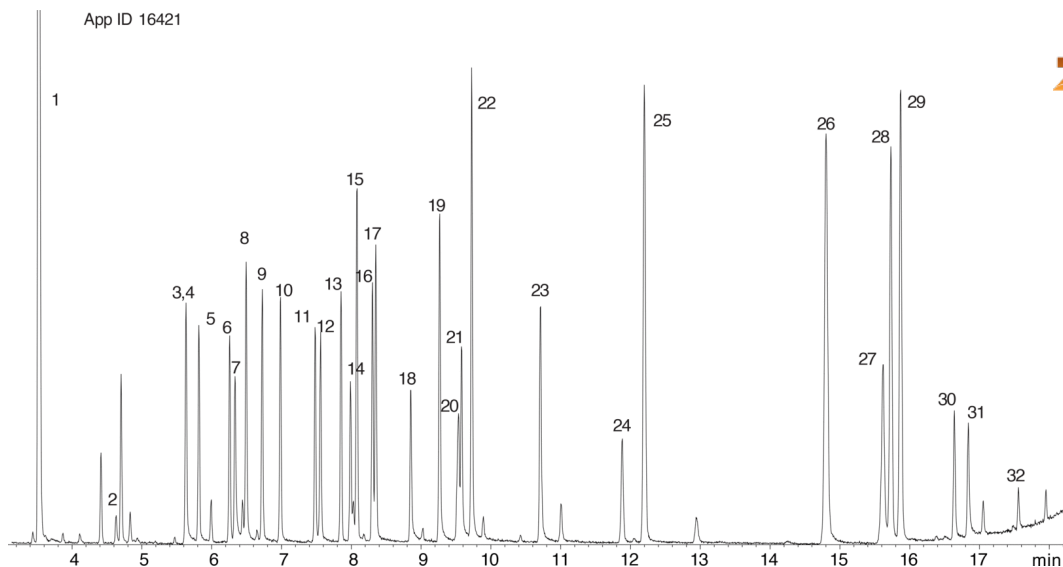


Drug screen on ZB-5plus

Column: Zebron ZB-5plus Capillary GC, Column 30m x 0.25mm x 0.25um, Ea
Phase: Proprietary Pesticides Phase
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G032-11
Oven Profile: 140 °C to 240 °C @ 10 °C/min for 5 min to 320 °C @ 25 °C/min for 2.25 min.
Carrier Gas: Constant Flow Helium, 1.1 mL/min
Injection: Split 15:1 1 µL @ 240°C
Detection: Mass Selective (MSD) (230°C)
Analyst Note: Analytes 25ppm in dichloromethane and are underivatized.

**Products used in this application:**

Drug screen on ZB-5plus

ANALYTES:

- 1 Nicotine
- 2 Methylecgonine
- 3 Ibuprofen
- 4 Allobarbitol
- 5 Aprobital
- 6 Butabital
- 7 Acetaminophen
- 8 Phenacetin
- 9 Amobarbital
- 10 Pentobarbital
- 11 Secobarbital
- 12 Meprobamate
- 13 Methyl benzilate
- 14 Caffeine
- 15 Benzphetamine
- 16 Hexobarbital
- 17 Dimenhydrinate
- 18 Doxylamine
- 19 Phenobarbital
- 20 8-Chlorotheophylline
- 21 Methapyrilene
- 22 Chlorpheniramine
- 23 Brompheniramine
- 24 Cocaine
- 25 Chlorcyclizine
- 26 Codeine
- 27 Morphine
- 28 Diazepam
- 29 Hydrocodone
- 30 6-Monoacetylmorphine
- 31 Oxycodone
- 32 Diacetylmorphine (Heroin)

