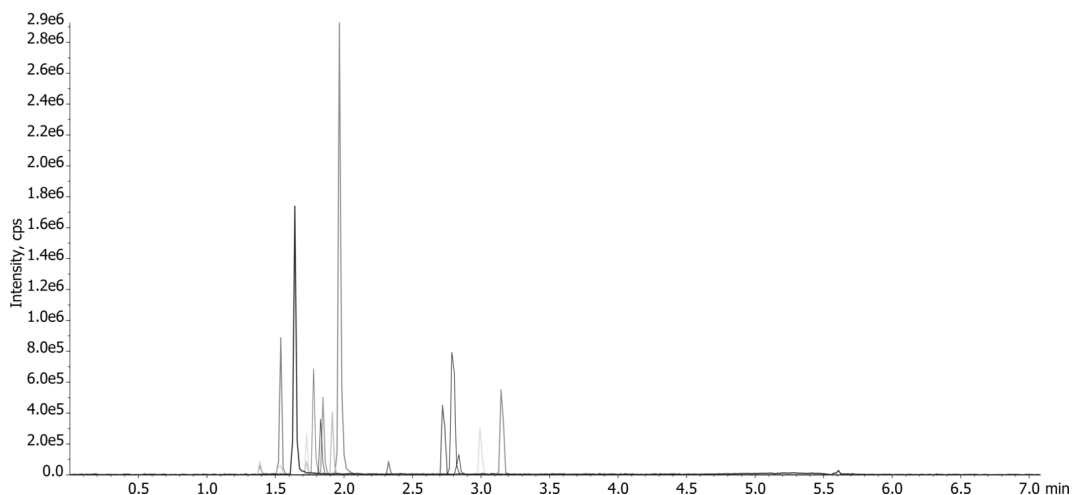


strata basic screen**Column:** Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 3.0 mm, Ea**Dimensions:** 50 x 3 mm ID**Order No:** 00B-4622-Y0**Elution Type:** Gradient**Eluent A:** 0.1% formic acid in water**Eluent B:** 0.1% formic acid in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	4	0	100
	3	5	0	100
	4	5.01	95	5
	5	7.01	95	5

Flow Rate: 600 µL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)

24639

**Products used in this application:**

strata basic screen

ANALYTES:

- 1** 6-MAM
- 2** Amphetamine
- 3** benzoylecgonine
- 4** buprenorphine
- 5** codeine
- 6** diazepam
- 7** fentanyl
- 8** hydrocodone
- 9** hydromorphone
- 10** methamphetamine
- 11** morphine
- 12** norbuprenorphine
- 13** nordiazepam
- 14** oxycodone
- 15** PCP
- 16** Temazepam
- 17** oxazepam

