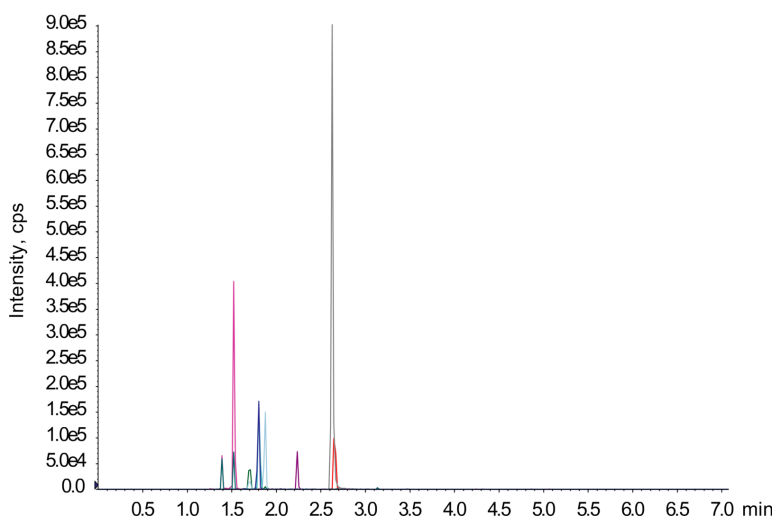


opiates on strata basic screen**Column:** Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 3.0 mm, Ea**Dimensions:** 50 x 3 mm ID**Order No:** 00B-4622-Y0**Elution Type:** Gradient**Eluent A:** 0.1% formic acid in water**Eluent B:** 0.1% formic acid in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	4	0	100
	3	4.5	0	100
	4	7.1	95	5

Flow Rate: 600 µL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)
app ID 24664**Products used in this application:****ANALYTES:****1** 6-MAM

Retention Time: 1.8 min

2 Buprenorphine

Retention Time: 2.66 min

3 codeine

Retention Time: 1.71 min

4 fentanyl

Retention Time: 2.63 min

5 hydrocodone

Retention Time: 1.88 min

6 hydromorphone

Retention Time: 1.53 min

7 morphine

Retention Time: 1.39 min

8 norbuprenorphine

Retention Time: 2.23 min

9 oxycodone

Retention Time: 1.82 min

