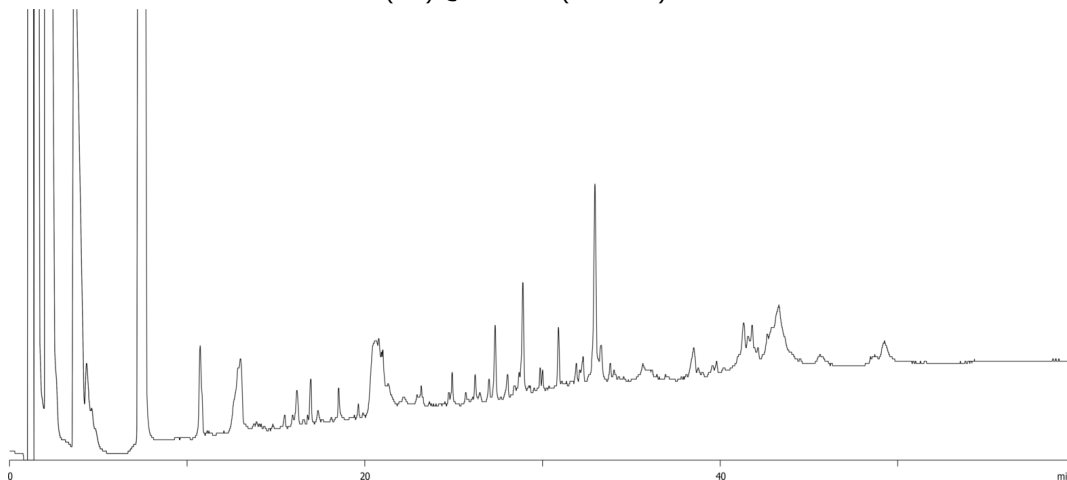


Anti Rituximab Fab Peptide Mapping, bioZen Peptide PS-C18**Column:** bioZenTM 1.6 μ m Peptide PS-C18, LC Column 150 x 2.1 mm, Ea**Dimensions:** 150 x 2.1 mm ID**Order No:** 00F-4770-AN**Elution Type:** Gradient**Eluent A:** 0.1%TFA in water**Eluent B:** 0.1% TFA in acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	99	1
	2	5	99	1
	3	60	50	50
	4	65	50	50
	5	70	5	95
	6	75	5	95
	7	78	5	95

Flow Rate: 0.3 mL/min**Col. Temp.:** 40 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)**ANALYTES:****1** anti Rituximab FAB