

Trastuzumab mcMMAF on bioZen 1.8 µm SEC-3

Column: bioZen™ 1.8 µm SEC-3 , LC Column 300 x 4.6 mm, Ea

Dimensions: 300 x 4.6 mm ID

Order No: 00H-4772-E0

Elution Type: Isocratic

Eluent A: 50 mM K₂PO₄ 100 mM Na₂SO₄, pH 6.8□□

Gradient


Step No.	Time (min)	Pct A
----------	------------	-------

Profile:	1	30	100
-----------------	----------	----	-----

Flow Rate: 0.2 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (25 °C)

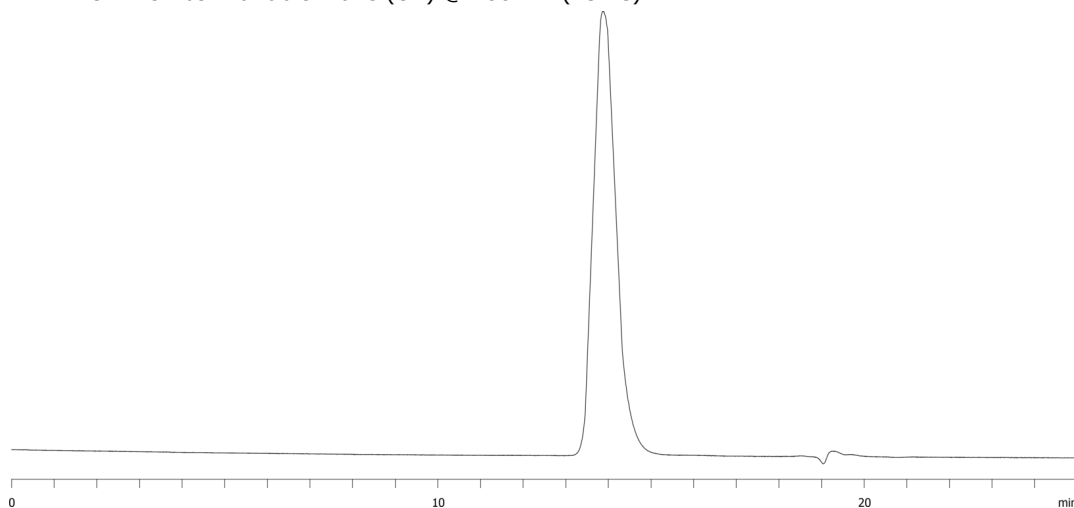

ANALYTES:

1 MMAE

