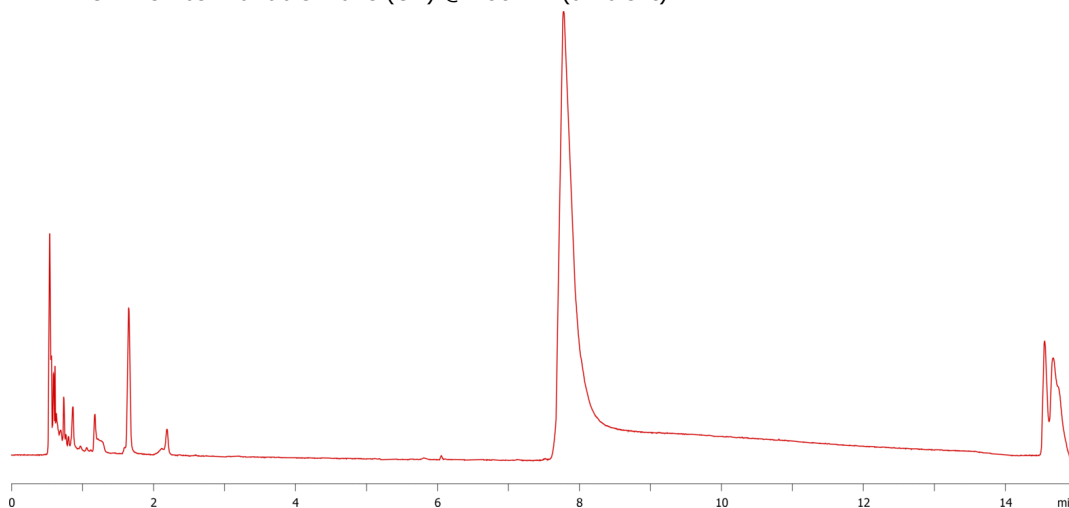


Trastuzumab Intact 80°C**Column:** bioZen[™] 3.6 µm Intact XB-C8, LC Column 150 x 2.1 mm, Ea**Dimensions:** 150 x 2.1 mm ID**Order No:** 00F-4766-AN**Elution Type:** Gradient**Eluent A:** 0.1% TFA in water**Eluent B:** 0.1% TFA in acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	80	20
	2	1	80	20
	3	13	55	45

Flow Rate: 0.5 mL/min**Col. Temp.:** 80 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)**ANALYTES:****1** Trastuzumab