

bioZen SEC-3 Batch To Batch Reproducibility

Column: bioZen™ 1.8 µm SEC-3 , LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4772-E0

Elution Type: Isocratic


Eluent A: 100 mM Sodium Phosphate, pH 6.8

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.3 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

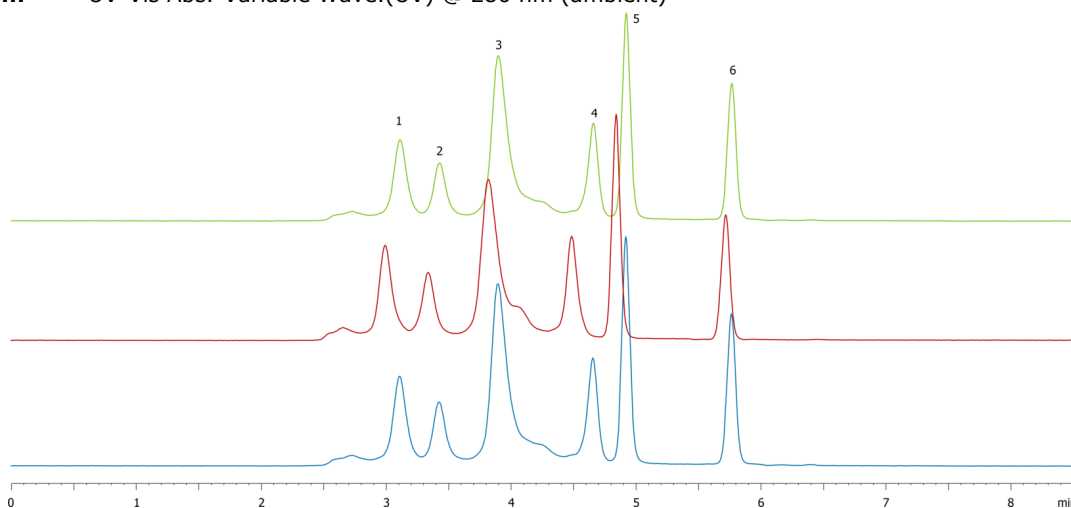

ANALYTES:

- 1 . Thyroglobulin (669 kDa)
- 2 IgA (300 kDa)
- 3 IgG (150 kDa)
- 4 . Ovalbumin (44 kDa)
- 5 Myoglobin (17 kDa)
- 6 Uridine

