

Zebon ZB-624PLUS QC Test Mix

Column: Zebon ZB-624plus, GC Cap. Column 30m x 0.32mm x 1.8µm, Ea

Phase: Proprietary Pesticides Phase

Dimensions: 30 meters x 0.32 mm x 1.8 µm

Order No: 7HM-G040-31

Oven Profile: 180°C for 12 min


Carrier Gas: Constant Flow Hydrogen, 40 cm/sec

Injection: Split 100:1 1 µL @ 250°C

Detection: GC/FID (305°C)



Products used in this application:

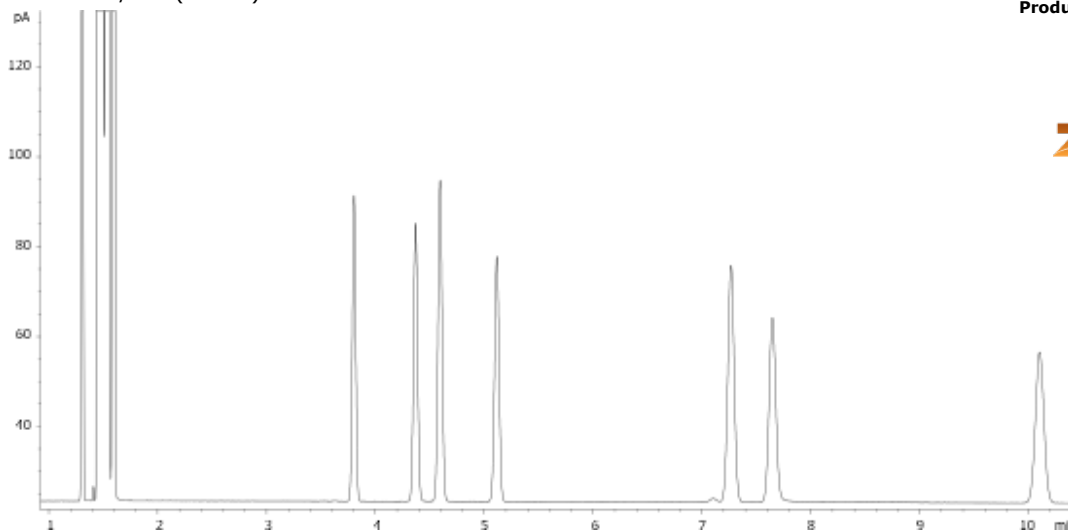

ANALYTES:

- 1 Dodecane
- 2 2,4-Dimethylphenol
- 3 2,6-Dimethylaniline
- 4 Tridecane
- 5 1-Methylnaphthalene
- 6 1-Undecanol
- 7 Pentadecane
- 8

