

16 Sugars on Luna Omega 3µm Sugar 250x4.6mm using ELSD

Column: Luna® Omega 3 µm SUGAR 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4775-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Acetonitrile


Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	70	20	80

Flow Rate: 2 mL/min

Col. Temp.: 25 °C

Detection: LC/ELSD (Light Scattering) @ (ambient)

Analyst Note: Water:Acetonitrile (20:80)

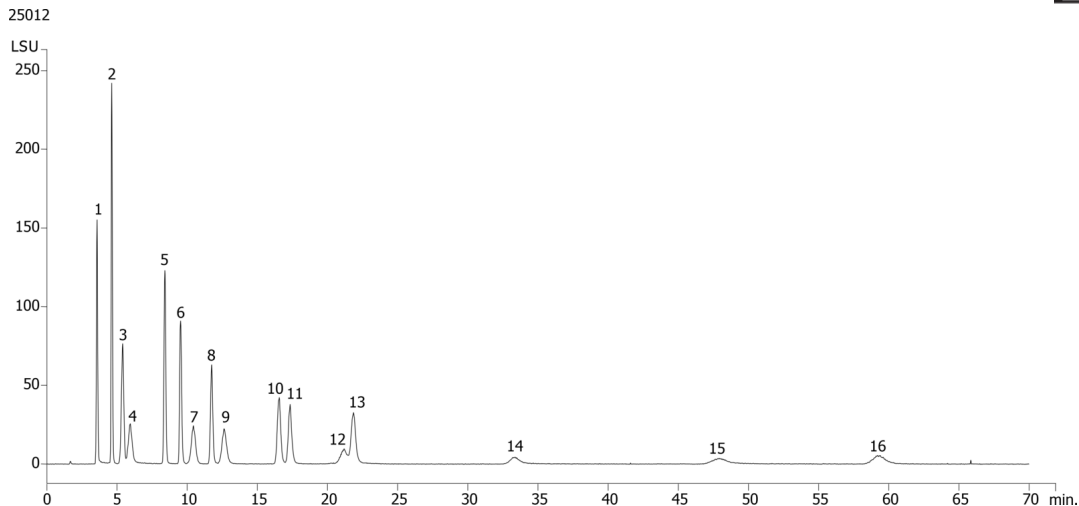

Products used in this application:



16 Sugars on Luna Omega 3µm Sugar 250x4.6mm using ELSD

ANALYTES:

- 1** Xylose
- 2** Fructose
- 3** Glucose
- 4** Galactose
- 5** Sucrose
- 6** Turanose
- 7** Maltose
- 8** Trehalose
- 9** Lactose
- 10** Erllose
- 11** Melezitose
- 12** Maltotriose
- 13** Raffinose
- 14** Isomaltotriose DP3
- 15** Cellotetraose
- 16** Stachyose

