

16 Cannabinoids on Kinetex C18

Column: Kinetex® 2.6 µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4462-E0

Elution Type: Gradient

Eluent A: 0.1% formic acid in water

Eluent B: 0.1% formic acid in 75:25 Methanol:Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	20	80
	2	8	20	80
	3	10	0	100
	4	12	0	100


Flow Rate: 1.5 mL/min

Col. Temp.: 30 °C

Detection: LC/UV (DAD, PDA) @ 228.0000000000 (ambient)



Products used in this application:

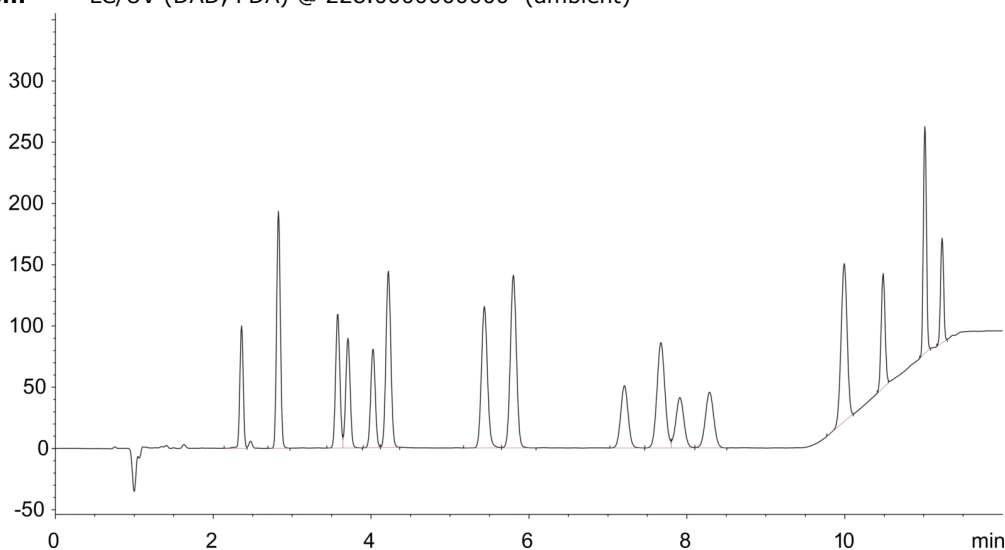

16 Cannabinoids on Kinetex C18

ANALYTES:

- 1** CBDV
- 2** CBDVA
- 3** CBD
- 4** CBG
- 5** THCV
- 6** CBDA
- 7** CBGA
- 8** CBN
- 9** D9-THC
- 10** THCVA
- 11** D8-THC
- 12** CBL
- 13** CBC
- 14** CBNA
- 15** THCA-A
- 16** CBCA

