

Lux i-Amylose-3 Separation for Thalidomide

Column: Lux® 5 µm i-Amylose-3, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4779-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Diethylamine

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	14	100	0.1

Flow Rate: 1 mL/min


Col. Temp.: 25 °C

Detection: LC/UV (DAD, PDA) @ 254.000000000 nm (nanometers) (ambient)

Analyst Note: 2 mg/mL in Ethanol



Products used in this application:

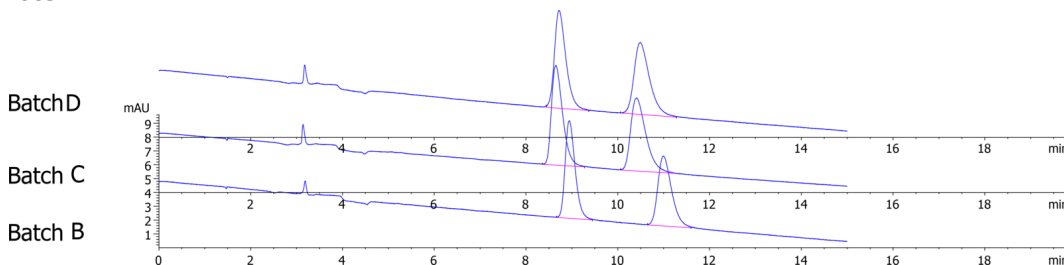

ANALYTES:

- 1 Thalidomide
Retention Time: 8.944 min
- 2 Thalidomide
Retention Time: 10.995 min

