

ASTM D5441 (High Concentration Mix-Sigma Aldrich 47942) on ZB-DHA-PONA 50m x 0.20mm x 0.50µm

Column: Zebron ZB-DHA-PONA, GC Cap. Column 50m x 0.20mm x 0.5µm , Ea

Phase:

Dimensions: 50 meters x 0.2 mm x 0.5 µm

Order No: 7JE-G042-17

Oven Profile: 40°C for 13 min to 180°C @ 10°C/min 180°C for 3 min

Carrier Gas: Constant Flow Helium, 0.9 mL/min


Injection: Split 50:1 0.5 µL @ 200°C

Detection: GC/FID (250°C)

Analyst Note: Resolution between trans-2-pentene/tert-butanol and tert-butanol/cis-2-pentene meets the ASTM Method D5441



Products used in this application:

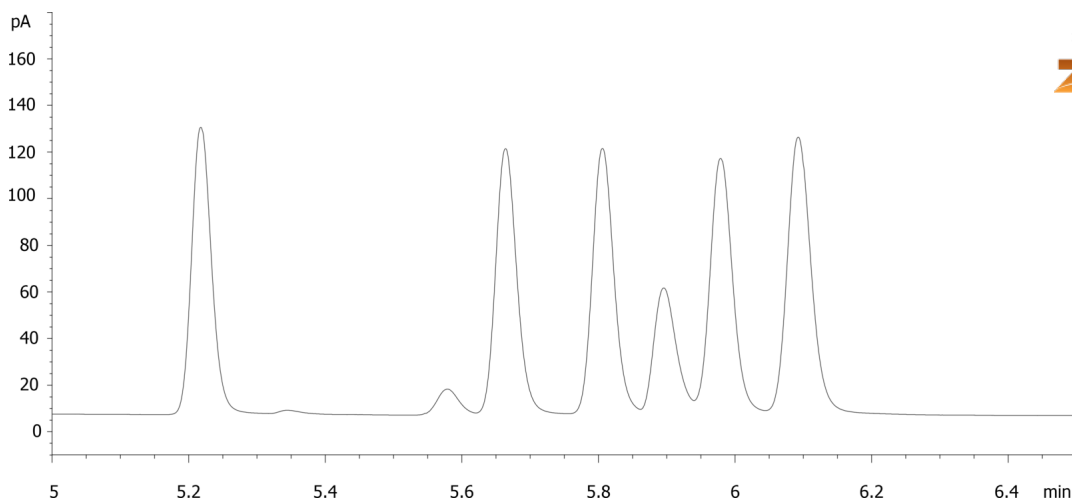

ASTM D5441 (High Concentration Mix-Sigma Aldrich 47942) on ZB-DHA-PONA 50m x 0.20mm x 0.50µm

ANALYTES:

- 1** Methanol
- 2** Iso-Pentane
- 3** n-Pentane
- 4** Trans-2-Pentene
- 5** Tert-Butanol
- 6** Cis-2-Pentene
- 7** 2-Methyl-2-Butene
- 8** MTBE
- 9** Tert-Butyl Ethyl ether
- 10** Tert-Amyl methyl ether
- 11** 2,4,4-Trimethyl-1-pentene
- 12** 4,4--Dimethyl-2-neopenyl-1-pentene (triisobutylene)
- 13** 2,2,4,6,6-Pentamethyl-3-heptene (triisobutylene)

