

Anticonvulsants

Column: Luna® 5 µm C8(2) 100 Å, LC Column 150 x 2 mm, Ea

Dimensions: 150 x 2 mm ID

Order No: 00F-4249-B0

Elution Type: Isocratic

Eluent A: 20 mM Ammonium acetate

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.2 mL/min

Col. Temp.: ambient

Detection: Mass Spectrometer (MS) @ nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 2.0mm, 10/Pk Part No.: AJ0-4289

- Holder Part No.: KJ0-4282



Products used in this application:

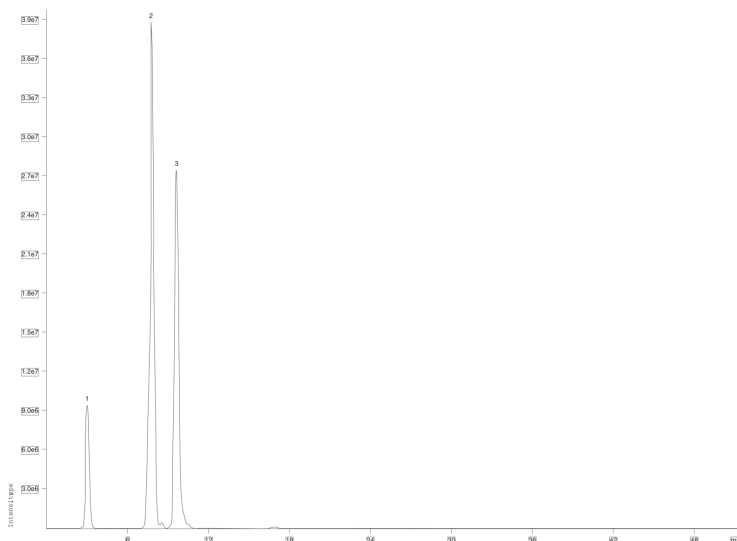

ANALYTES:

- 1 Primidone
- 2 Carbamazepine
- 3 Phenytoin

