

Ivermectin and Related Substances Test Solution - Ph. Eur. Monograph 1336

Column: Luna® 5 µm C18(2) 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4252-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Methanol

Eluent C: Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B	Pct C
----------	----------	------------	-------	-------	-------

Profile:	1	0	15	34	51
-----------------	----------	---	----	----	----

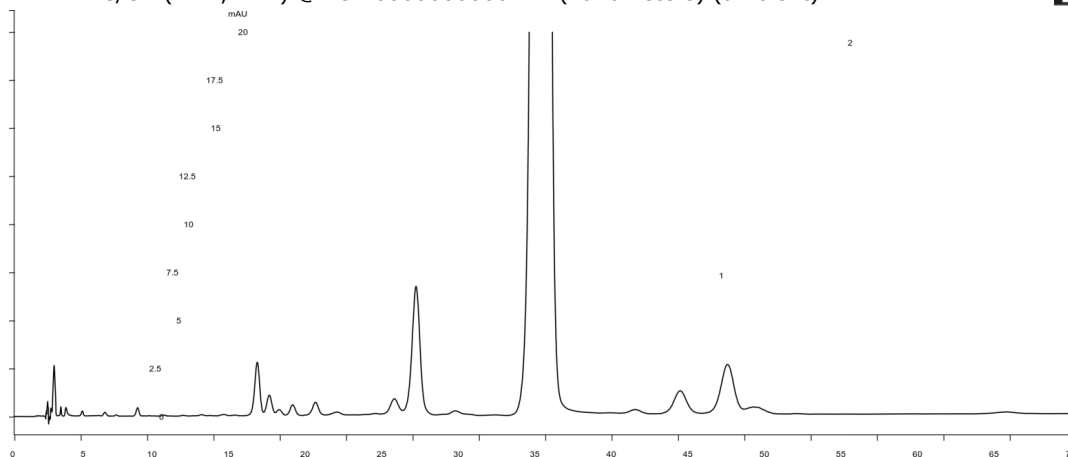
Flow Rate: 1 mL/min

Col. Temp.: 25 °C

Detection: LC/UV (DAD, PDA) @ 254.0000000000 nm (nanometers) (ambient)



Products used in this application:

**ANALYTES:**

1 H2B1b

Retention Time: 30.275 min

2 H2B1a

Retention Time: 39.57 min

