

Rituximab on bioZen 2.6 µm WidePore C4

Column: bioZen™ 2.6 µm WidePore C4, LC Column 100 x 2.1 mm, Ea

Dimensions: 100 x 2.1 mm ID

Order No: 00D-4786-AN

Elution Type: Gradient

Eluent A: 0.1 % TFA in Water


Eluent B: 0.1% TFA in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30
	2	5	60	40
	3	6	20	80
	4	6.1	70	30
	5	10	70	30

Flow Rate: 0.8 mL/min

Col. Temp.: 60 °C

Detection: LC/UV (DAD, PDA) @ (ambient)

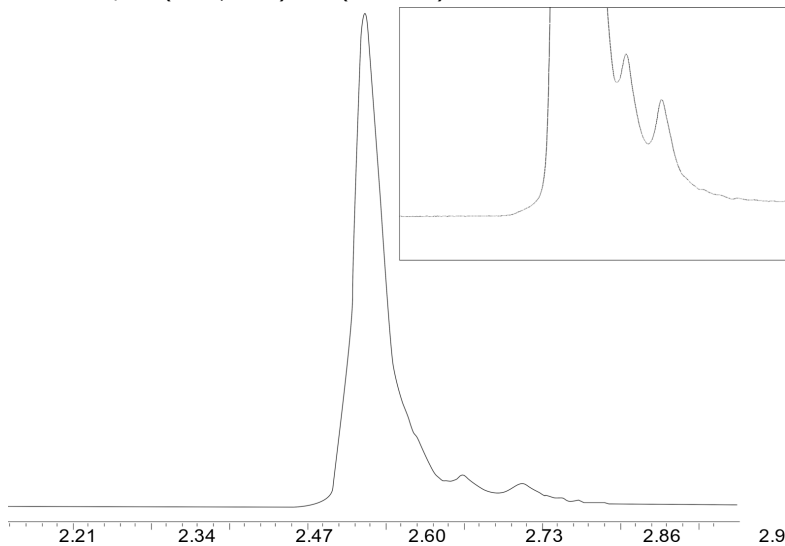

ANALYTES:

