


TCDF Separation using Zebron ZB-Dioxin GC Column**Column:** Zebron™ ZB-Dioxin, GC Cap. Column 60m, 0.25mm, 0.20um, Ea**Phase:****Dimensions:** 60 meters x 0.25 mm x 0.2 µm**Order No:** 7KG-G045-10**Oven Profile:** 160°C for 2.4 min to 200°C @ 25°C/min to 220°C @ 5°C/min 220°C for 19 min to 288°C @ 4°C/min to 300°C @ 5°C/min 300°C for 7.6 min**Carrier Gas:** Constant Flow Helium, 1.25 mL/min**Injection:** Splitless :1 1 µL @ 280°C**Detection:** GC/HRMS (300°C)**Products used in this application:****ANALYTES:**

- 1** 1,3,6,8-TCDF
Retention Time: 23.2 min
- 2** 1,3,4,7-TCDF
Retention Time: 29.57 min
- 3** 2,3,7,8-TCDF
Retention Time: 29.98 min
- 4** 1,2,3,9-TCDF
Retention Time: 31.14 min
- 5** 1,2,8,9-TCDF
Retention Time: 32.79 min

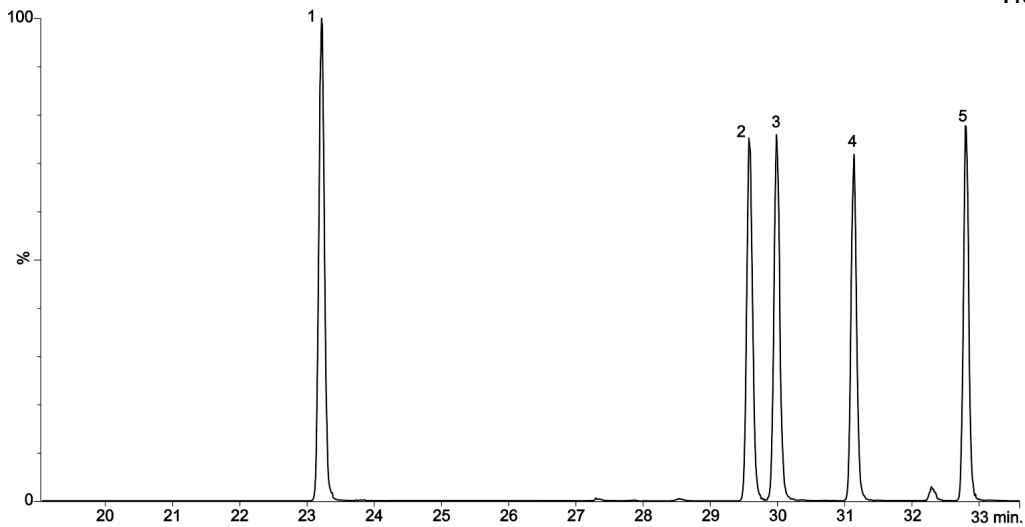