

HPLC of Carbohydrates from Yellow-Seeded Canola on Rezex RSO

Column: Rezex™ RSO-Oligosaccharide Ag+ 4%, LC Column 200 x 10 mm, Ea

Dimensions: 200 x 10 mm ID

Order No: 00P-0133-N0

Elution Type: Isocratic

Eluent A: Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient


Detection: Refractive Index (RI) @ nm (60 °C)
2696

REZEX™
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX™

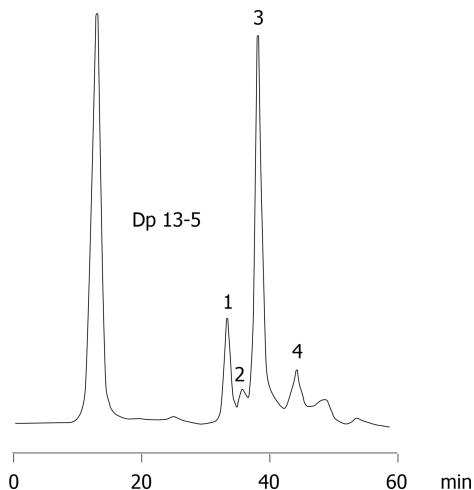

ANALYTES:

- 1 Stachyose
- 2 Raffinose
- 3 Sucrose
- 4 Glucose

