

Lorazepam Tablets - USP Method

Column: Luna® 5 µm C18(2) 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4252-E0

Elution Type: Isocratic

Eluent A: Water

Eluent B: Methanol

Eluent C: Acetic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
1	0	54	44	2	

Flow Rate: 2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282



Products used in this application:

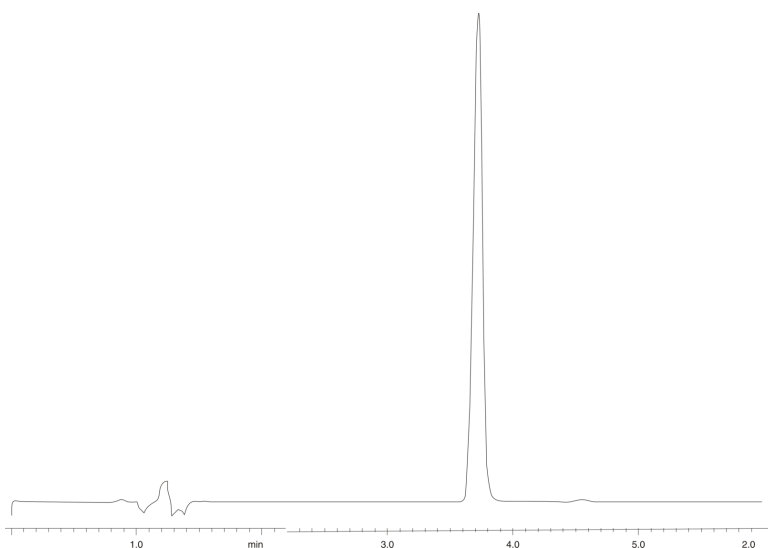

ANALYTES:

1 Lorazepam

Retention Time: 3.709 min

